


# Vascular Technology

AN ILLUSTRATED REVIEW

5th EDITION

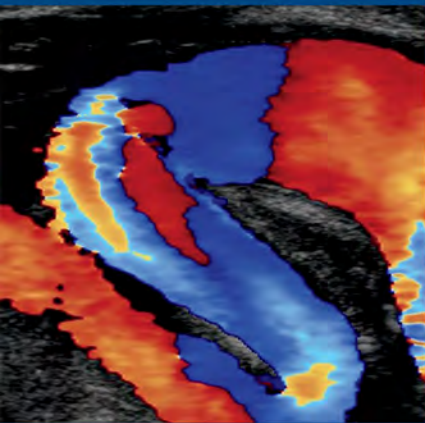

3



## 1-2-3 Step Ultrasound Education & Test Preparation

### Step 1


Review text

### Step 2

Mock examination

### Step 3

Q&A memory skills  
flashcard drill

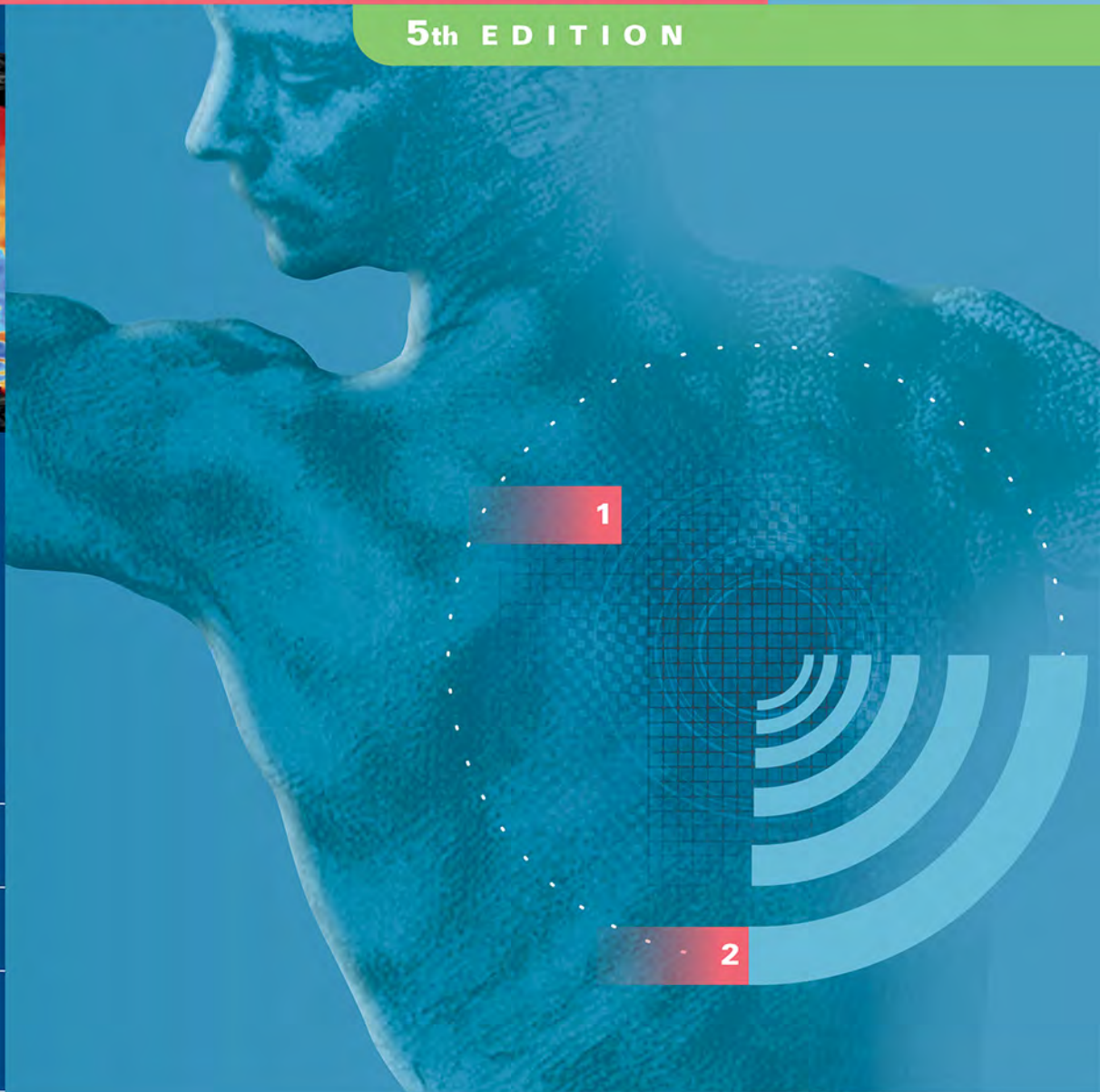

1

2

SDMS-Approved  
Continuing Education Activity

Approved for **15** Hours CME Credit

CLAUDIA RUMWELL

MICHAELNE McPHARLIN

# Vascular Technology

An Illustrated Review

5<sup>th</sup> Edition

**Claudia Rumwell**, RN, RVT, FSVU

**Michalene McPharlin**, RN, RVT, RVS, FSVU





Copyright © 2015 by Davies Publishing, Inc.

All rights reserved. No part of this work may be reproduced, stored in a retrieval system, or transmitted in any form or by any means, electronic or mechanical, including photocopying, scanning, and recording, without prior written permission from the publisher.



**Davies Publishing, Inc.**

*Specialists in Ultrasound Education,  
Test Preparation, and Continuing  
Medical Education*

32 South Raymond Avenue  
Pasadena, California 91105-1961

Phone 626-792-3046

Facsimile 626-792-5308

e-mail [info@daviespublishing.com](mailto:info@daviespublishing.com)

[www.daviespublishing.com](http://www.daviespublishing.com)

Michael Davies, Publisher

Charlene Locke, Production Manager

Christina J. Moose, Editorial Director

Janet Heard, Operations Manager

Daniel Liota, Digital Media

Satori Design Group, Inc., Design

Denise Eggman, Illustration

Stephen Beebe, Illustration

*Notice to Users of This Publication:*

In the field of ultrasonography, knowledge, technique, and best practices are continually evolving. With new research and developing technologies, changes in methodologies, professional practices, and medical treatment may become necessary.

Sonography practitioners and other medical professionals and researchers must rely on their experience and knowledge when evaluating and using information, methods, technologies, experiments, and medications described herein, always remaining mindful of their own, their patients', their coworkers', and others' safety and well-being.

Regarding any treatments, procedures, technologies, and/or pharmaceutical products identified, users of this publication are advised to check the most current information provided by product manufacturers and professional societies to verify their latest recommendations regarding methodologies, dosages, timing of administration, practice guidelines, standards of care, and contraindications. It is the responsibility of practitioners, relying on their own experience and knowledge of their patients, to report, make recommendations, and work within these standards, in consultation with other clinicians such as referring and supervisory physicians, to determine to take all appropriate safety procedures and to determine the best treatment of each individual patient.

To the fullest extent of the law, neither Davies Publishing, Inc., nor the authors, contributors, reviewers, or editors assume any liability for any injury and/or damage to persons or property as a matter of products liability, negligence, or otherwise, or from any use or operation of any methods, products, instructions, or ideas contained in the material herein.

**Library of Congress Cataloging-in-Publication Data**

Rumwell, Claudia, 1946- author.

Vascular technology : an illustrated review / Claudia Rumwell, Michalene McPharlin. — 5th edition.  
p. ; cm.

Includes bibliographical references and index.

ISBN 978-0-941022-85-9 (alk. paper) — ISBN 0-941022-85-4 (alk. paper)

I. McPharlin, Michalene, 1949- author. II. Title.

[DNLM: 1. Vascular Diseases—diagnosis—Examination Questions. 2. Diagnostic Techniques, Cardiovascular—Examination Questions. WG 18.2]

RC691.6.A53

616.1'30076—dc23

2014014684

Printed and bound in China

ISBN 978-0-941022-85-9

# Contents

Reviewers vii

Preface to the 5th Edition ix

Acknowledgments xi

## **PART I ARTERIAL EVALUATION 1**

**Chapter 1** Gross Anatomy of the Central and Peripheral Arterial Systems 3

**Chapter 2** Physiology and Fluid Dynamics 15

**Chapter 3** Testing Considerations, Patient History, Mechanisms of Disease, and Physical Examination 43

**Chapter 4** Doppler Waveform Analysis in the Upper and Lower Extremities 57

**Chapter 5** Doppler Segmental Pressures—Lower Extremities 69

**Chapter 6** Doppler Segmental Pressures—Upper Extremities 81

**Chapter 7** Laser Doppler 89

**Chapter 8** Penile Pressures and Waveforms—Duplex/Color Flow Imaging Evaluation 93

**Chapter 9** Plethysmography—Upper and Lower Extremities 99

**Chapter 10** Digital Pressures and Plethysmography 107

**Chapter 11** Transcutaneous Oximetry (TcPO<sub>2</sub>) 115

**Chapter 12** Duplex Scanning and Color Flow Imaging of the Upper Extremities 119

**Chapter 13** Duplex Scanning and Color Flow Imaging of the Lower Extremities 127

**Chapter 14** Duplex Scanning and Color Flow Imaging of the Abdominal Vessels 139

**Chapter 15** Preoperative Mapping Procedures 159

**Chapter 16** Atypical Vascular Disorders 165

**Chapter 17** Alternative Diagnostic Tests and Therapeutic Interventions 173

**PART II CEREBROVASCULAR EVALUATION 197**

**Chapter 18** Gross Anatomy, Physiology, and Fluid Dynamics 199

**Chapter 19** Testing Considerations, Patient History, Mechanisms of Disease, and Physical Examination 207

**Chapter 20** Carotid Duplex Scanning and Color Flow Imaging 221

**Chapter 21** Transcranial Doppler (TCD) 239

**Chapter 22** Atypical Vascular Disorders 249

**Chapter 23** Alternative Diagnostic Tests and Therapeutic Interventions 253

**PART III VENOUS EVALUATION 263**

**Chapter 24** Gross Anatomy of the Central and Peripheral Venous Systems 265

**Chapter 25** Venous Hemodynamics 277

**Chapter 26** Testing Considerations, Patient History, Mechanisms of Disease, and Physical Examination 281

**Chapter 27** Venous Photoplethysmography 293

<b>Chapter 28</b>	Venous Air Plethysmography	297
<b>Chapter 29</b>	Continuous-Wave Doppler	301
<b>Chapter 30</b>	Duplex Scanning and Color Flow Imaging in Venous Evaluation	307
<b>Chapter 31</b>	Alternative Diagnostic Tests	331
<b>Chapter 32</b>	Medical Therapies	337
<b>Chapter 33</b>	Surgical and Endovascular Therapies and Other Treatment Options	341
<b>PART IV</b>	<b>TEST VALIDATION, STATISTICS, AND PATIENT SAFETY</b>	<b>353</b>
<b>Chapter 34</b>	Statistical Profile and Correlation	355
<b>Chapter 35</b>	Patient Care and Environmental Issues	363
Appendix A	Case Studies for Self-Assessment	367
Appendix B	Application for CME Credit	401
Appendix C	Bibliography	433
Index		435

# Physiology and Fluid Dynamics

*Arterial System*

*Energy*

*Blood Flow Characteristics*

*Poiseuille's Law*

*Reynolds Number*

*Pressure/Velocity Relationships (Bernoulli Principle)*

*Steady versus Pulsatile Flow*

*Peripheral Resistance*

*Cardiac Effects*

*Collateral Effects*

*Effects of Exercise*

*Effects of Stenosis on Flow*

## ARTERIAL SYSTEM


The arterial system is a multibranched elastic conduit that carries blood away from the heart and outward to the most distant tissues.

- The arterial tree oscillates with every beat of the heart, each one of which pumps approximately 70 ml of blood into the aorta and causes a blood pressure pulse.
- At the beginning of the cardiac contraction, the pressure in the left ventricle rises rapidly, quickly exceeding that in the aorta so that the aortic valve opens, blood is ejected, and the blood pressure rises. The amount of blood ejected is called the *stroke volume*.
- Increased heart rate delivers an increased blood volume that supplies more nutrients. Conversely, the lower the heart rate, the smaller the volume of pumped blood.
- The patient's cardiac status plays an important role in the movement of blood throughout the vascular system.

The heart pump generates the pressure (*potential energy*) to move the blood. The stroke volume of blood produced by each heartbeat creates a pressure (or *energy*) wave that travels rapidly throughout the arterial system (Figure 2-1):

- The propagation speed, shape, and strength of the pressure wave change as the wave moves through the arterial system.

**Figure 2-1.** Once the heart generates the pressure to move the blood, the energy wave produced travels rapidly throughout the system beginning with the left ventricle (A) to the capillary bed (E) back to the heart via the venous system (F–I).

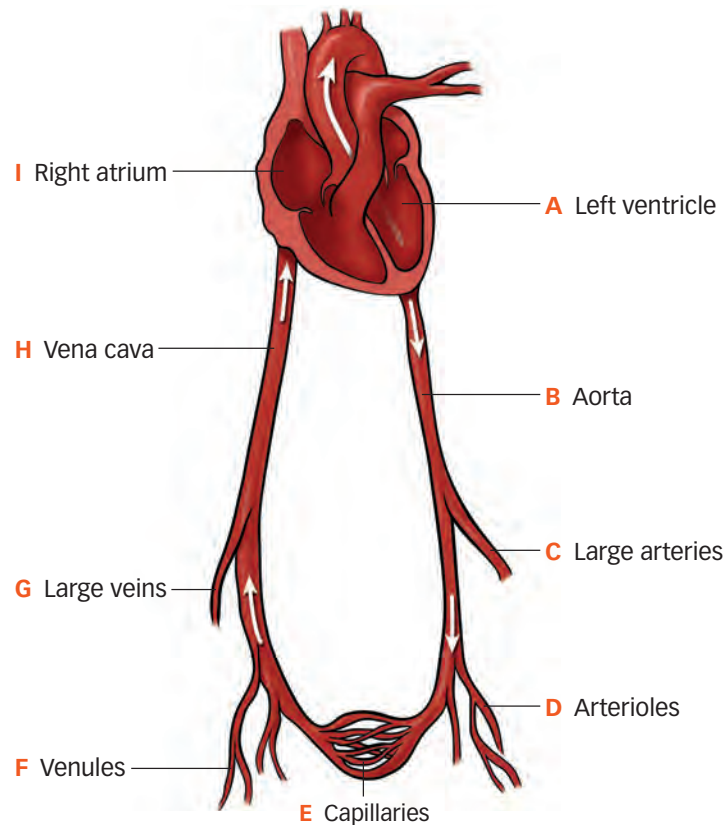



**Example:** As the arterial pressure wave moves distally, away from the heart and out toward the periphery, the propagation speed—the pulse wave velocity—increases with the growing stiffness of the arterial walls.

- Variations in the characteristics of the vessels influence these alterations in blood flow. Velocity and flow direction also vary with each heartbeat.
- As the pressure wave moves from the large arteries through the high-resistance vessels, capillaries, and then into the venous side, the mean pressure gradually declines because of losses in total fluid energy.

The pumping action of the heart maintains a high volume of blood in the arterial side of the system that in turn sustains a high pressure gradient between the arterial and venous sides of the circulation. This pressure gradient is necessary to maintain flow.

Cardiac output governs the amount of blood that enters the arterial system, while arterial pressure and total peripheral resistance (which is controlled by the level of vasoconstriction in the microcirculation) determine the volume of blood that leaves it.

A large portion of the energy created with each left ventricular contraction results in distention of the arteries, producing an arterial “reservoir” that stores some of the blood volume and the potential energy supplied to the system.

It is this store of energy and volume that promotes the flow of blood into the tissue during diastole. That is, potential energy is stored in the distended arterial wall and is released when the wall recoils.

Pressure is greatest at the heart and gradually decreases as the blood moves distally. This pressure difference (or *gradient*) is necessary to maintain blood flow.

## **ENERGY**


The movement of any fluid medium between two points requires two things: (1) a route along which the fluid can flow and (2) a difference in energy (pressure) levels between the two points. The volume of flow depends on the net energy difference between these two points, a factor that is affected by losses resulting from the movement of the fluid—i.e., friction—and any resistance within the pathway that opposes such movement.

The greater the energy difference (or the lower the resistance), the greater the flow, as illustrated in Figure 2-2.

The greater the pathway resistance and/or energy losses, the lower the flow.

***lower resistance = higher flow rate***  
***higher resistance = lower flow rate***


**Figure 2-2.** The greater the pressure gradient, the greater the flow. Example **A** has less pressure difference and therefore less flow than example **B**, which has a larger pressure gradient and greater flow.

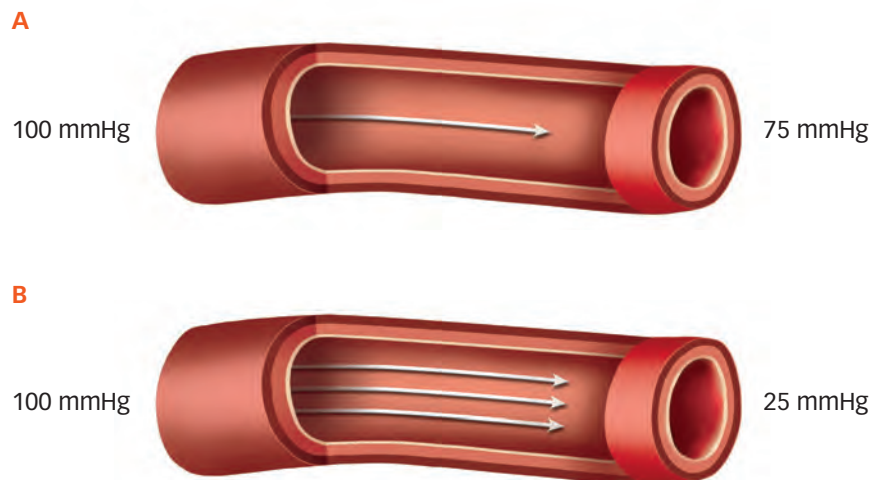

The total energy contained in moving fluid is the sum of pressure (potential) energy, kinetic (movement) energy, and gravitational energy, as explained below:

- *Pressure energy* is the main form of energy present in flowing blood (see Figure 2-3A). It is created by the pumping action of the heart, which subsequently distends the arterial vessels. This distention occurs because of elastin, which allows for a rapid increase in capacity. Arteries are more elastic proximally and less distally as they become smaller. Elastin quickly converts kinetic energy to pressure energy. The pressure energy is then converted back to kinetic energy in diastole. Pressure energy is also referred to as *potential energy*. Potential energy has several components. The dominant source is the pumping action of the heart muscle, as well as distention of the arterial wall.
- *Kinetic energy* is the ability of flowing blood to do work as a result of its velocity; it is the energy of something in motion (see Figure 2-3B). With regard to blood flow, the kinetic energy portion is small compared to the pressure energy. Kinetic energy is also proportional to the density of blood (which is normally stable) and to the square of its velocity. The everyday example of a dam illustrates the difference between potential and kinetic energy: The water behind the dam has *potential energy* (with the height of the dam providing a form of *gravitational energy*), while the water flowing through the dam has *kinetic energy*.
- *Gravitational energy (hydrostatic pressure)* is expressed in millimeters of mercury (mmHg). Changes in the height of the fluid column introduce the element of gravitational energy, which is hydrostatic pressure. In the circulatory system, hydrostatic pressure is equivalent to the weight of the column of blood extending from the heart, where the right atrium is considered the 0 pressure reference point (i.e., atmospheric pressure), to the level where the pressure is being measured.

**Figure 2-7. A** In a popliteal artery aneurysm, resistance to flow is reduced due to the larger-than-normal diameter. Although the volume of blood flow through this vessel remains constant, velocities are decreased compared to normal. **B** In this superficial femoral artery (SFA) stenosis, resistance to flow is increased because of the stenotic lumen. Although the volume flow remains constant, velocities are elevated as a result of the stenosis.

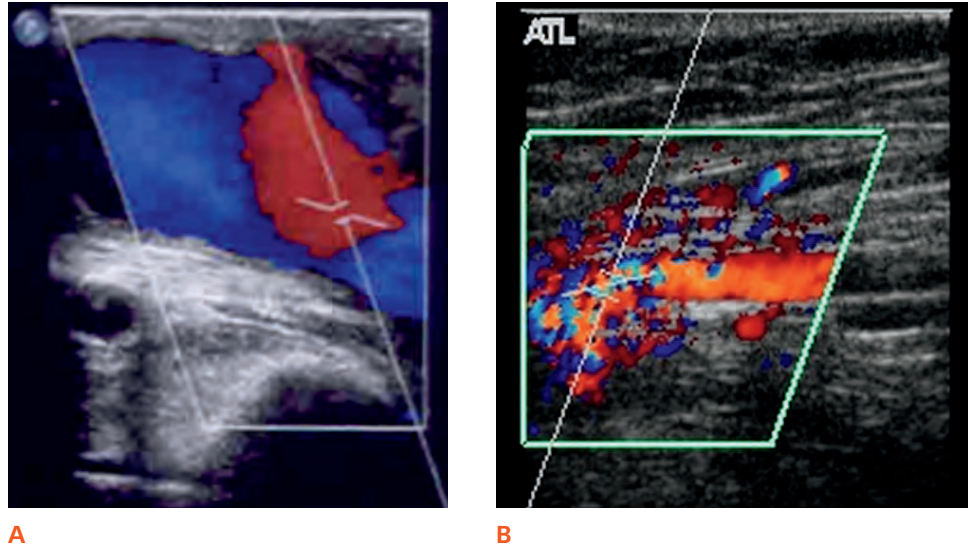

### POISEUILLE'S LAW


Poiseuille's law defines the relationship between volume flow (Q), pressure (P), and resistance (R) and may be written as  $Q = P/R$ . It helps to answer the question of how much fluid is moving through a vessel. When combined with the resistance equation (see page 21), Poiseuille's law may be stated as follows:

$$Q = \frac{(P_1 - P_2) \pi r^4}{8\eta L}$$

where Q = volume flow,  $P_1$  = pressure at the proximal end of the vessel,  $P_2$  = pressure at the distal end of the vessel, r = radius of the vessel, L = length of the vessel,  $\pi = 3.1416$ , and  $\eta$  = viscosity of the fluid.

- A change in the diameter (hence, radius) of a vessel affects resistance more dramatically than viscosity or vessel length (see Figures 2-7A and 2-7B).
- The radius of a vessel is directly proportional to the volume flow. Small changes in radius may result in large changes in flow, as demonstrated in Figure 2-8.

**Figure 2-8.** Parts **A** and **B** both illustrate narrowed segments. The reduction in volume flow through segment A will be greater than that through segment B, where the radius of the stenotic segment is somewhat larger.

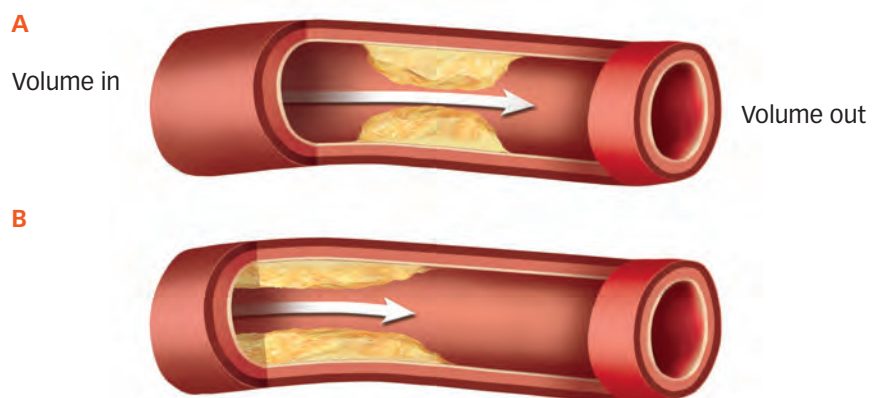

## REYNOLDS NUMBER

After the initial acceleration in systole, blood movement continues and develops into distinct streamline formations. Where  $P$  = pressure, note in Figure 2-10A that the streamlines are evenly distributed. When the flow pattern becomes unstable, these continuous streamlines break up and form small circular currents called eddy currents and vortices (swirling patterns of rotational flow) (Figure 2-10B).

Osborne Reynolds sought to determine how viscosity, vessel radius, and the pressure/volume relationship influence the stability of flow through a vessel. Although most of his work applied to straight, rigid tubes, it still provides insight into the physics of blood flow.

Flow volume increases as pressure increases, but only to a point. As flow changes from stable to disturbed, Reynolds found that an increase in pressure no longer increased flow volume. Instead, it increased flow disturbance, contributing to the formation of eddy currents.


The elements that affect the development of turbulent flow are expressed by a "dimensionless" number called the *Reynolds number* ( $Re$ ). The factors that affect the development of turbulence are expressed by this number according to the following equation:

$$Re = \frac{V\rho 2r}{\eta}$$

where  $Re$  = Reynolds number,  $V$  = velocity,  $\rho$  = the density of the fluid,  $r$  = the radius of the tube, and  $\eta$  = the viscosity of the fluid. Because the density and viscosity of the blood are fairly constant, the development of turbulence depends mainly on the size of the vessels and the velocity of flow. When the Reynolds number exceeds 2000, laminar flow tends to become disturbed. Flow disturbances also can occur at lower values because of other factors, such as body movement, pulsatility of blood flow, and irregularities of the vessel wall and plaque.

Turbulent flow may cause vessel walls to vibrate. The harmonics of this vibration produce vascular bruits.

**Figure 2-10.** Arterial flow streamlines. **A** Evenly distributed. **B** Disrupted (turbulent).  $P$  = pressure.

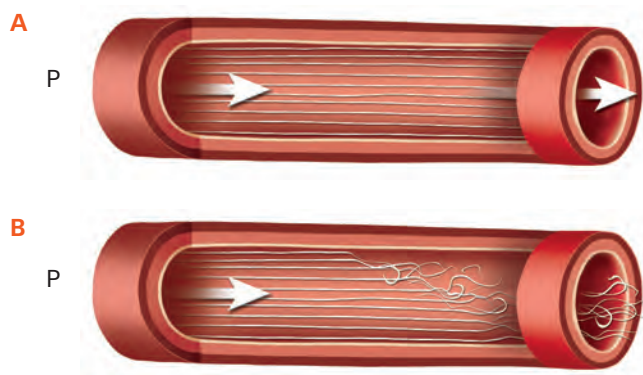

## PRESSURE/VELOCITY RELATIONSHIPS (BERNOULLI PRINCIPLE)

As previously described in the section on energy, the total energy contained in moving fluid is the sum of potential (i.e., pressure), kinetic, and gravitational energies. If one of these variables changes, the others also must change to maintain total fluid energy at the same level.


**Example:** If gravitational energy remains unchanged (that is, there is no change in the height of the fluid) but kinetic energy (velocity) increases, then potential (pressure) energy must decrease to maintain the same total fluid energy.

The Bernoulli equation shows that velocity and pressure are inversely related. Where there is high velocity, there is low pressure; where there is low velocity, there is high pressure. This inverse relationship between pressure and velocity explains why pressure decreases where fluid velocity increases (within the stenotic segment of an artery, for example) and why pressure distal to a stenosis (the region of poststenotic turbulence, where velocity decreases) is higher than that within the lesion itself. (See Figure 2-11.)

In other words, this is the law of conservation of energy:



1. In the region proximal to the stenosis (prestenosis), the pressure energy is higher and the kinetic energy lower. This region has the highest total energy sum.
2. As blood flows into the area of the stenosis, the pressure energy decreases and the kinetic energy increases (higher). However, total energy in this stenotic segment is less than that in the prestenotic segment because energy is lost (i.e., converted into heat) as it moves through the narrowing.
3. Lastly, distal to the stenosis (poststenosis), the kinetic energy decreases and the pressure energy increases. At this point the total energy sum is lowest.

Pressure gradients—the difference in pressure between two points in a vessel—are described as *flow separations*. Flow separations within a vessel may be

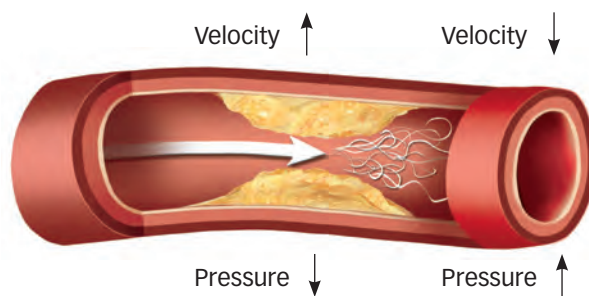

**Figure 2-11.** Velocity energy is elevated and pressure energy diminished within a stenosis. Relative velocity energy decreases and pressure energy increases distal to a stenosis. However, overall energy decreases.

**Figure 2-12.** Flow separation patterns. In both examples **A** and **B**, pressure energy is higher and velocity energy lower in the area of the pressure gradient (flow separation). This causes flow direction to move to the area of lower pressure energy. **A** Flow separation in the carotid bulb. **B** Flow separation at a curve (note that it is located on the inside of the vessel curve).

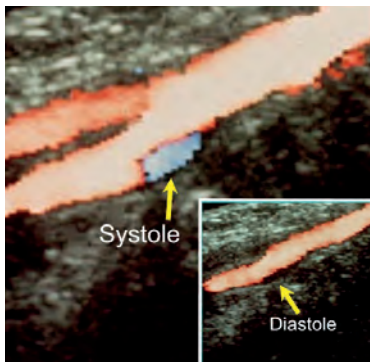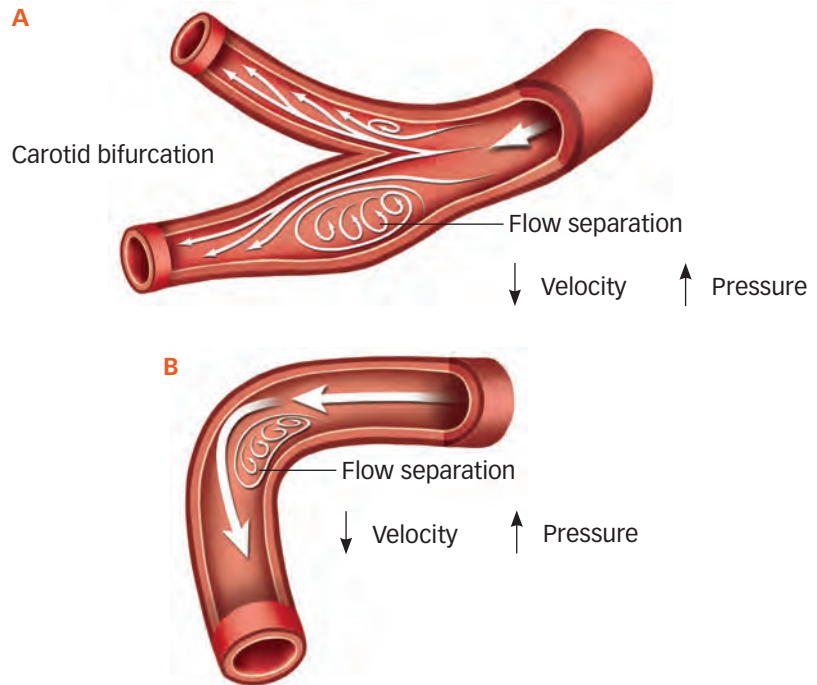

**Figure 2-13.** Longitudinal image of the carotid bifurcation in systole. The inset shows the same image in diastole. Note that flow separation (blue) is present during systole, absent in diastole.


caused by changes in the geometry of the vessel (with or without intraluminal disease) or the direction of the vessel, as depicted in Figures 2-12A and B. See also Figure 2-13.


Flow separations leave behind regions of flow reversal, stagnant or little movement. Figure 2-14 demonstrates two different regions within a bypass graft where these flow separations can occur.

Because flow moves from high to low pressure (described as a *pressure gradient*), the direction of flow in the region of flow separation (e.g., carotid bulb, bypass graft anastomosis) changes with respect to the transducer, causing a visible color change in the color flow image at systole.

During diastole, when flow at the vessel wall is stagnant, there is no movement of blood and therefore no color in the color flow image. The flow separation pattern is an ideal one to use to help define whether an image is in systole or diastole.

**Figure 2-14.** In this reversed saphenous vein graft at the end-to-side proximal anastomosis there is an area of flow stagnation (flow separation) on the inside wall. Distally, there are areas of flow separation at the site of a valve cusp. Note the variability of vessel dimension in the bypass graft. CFA = common femoral artery, SFA = superficial femoral artery.

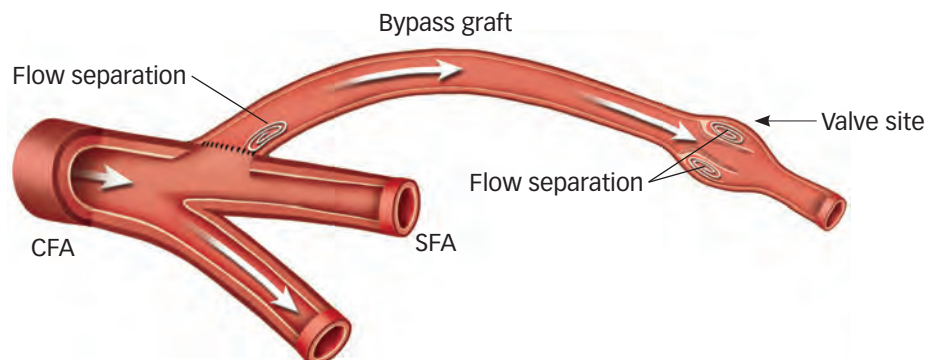

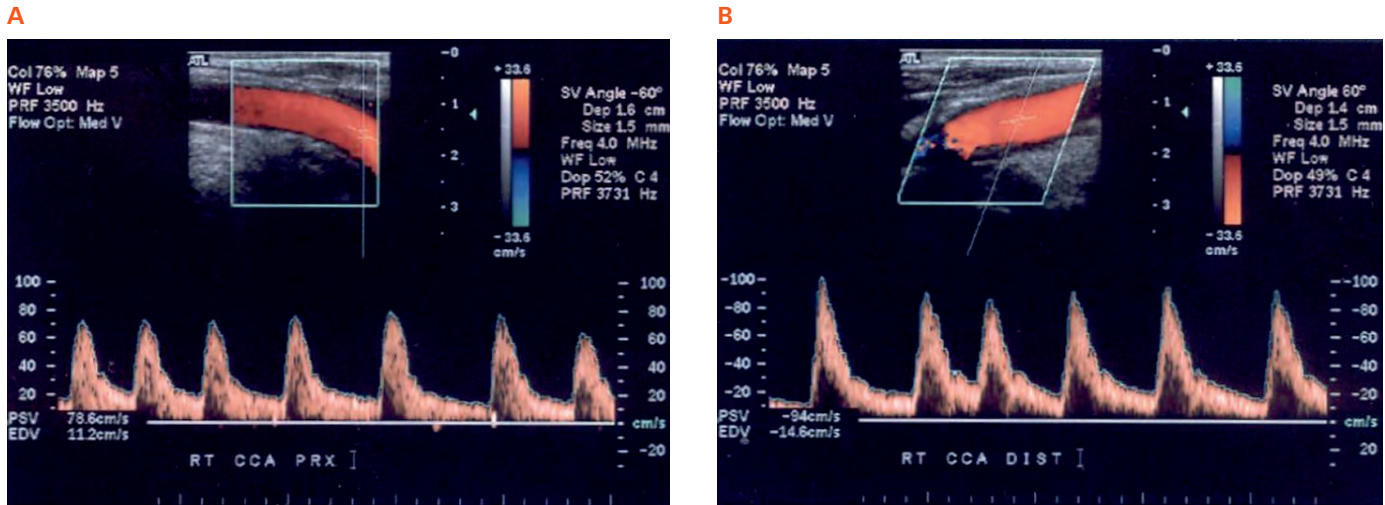


**Figure 2-22.** An irregular heart rhythm is evident in this spectral analysis of a common carotid artery.


- A severely irregular heart rhythm (Figure 2-23), on the other hand, represents a challenge in obtaining reliable peak systolic velocity (PSV) measurements. Some protocols call for averaging a few cycles to be the most reliable; others consider the averaging of 10 cycles to be the best method. When determining the significance of a stenosis, calculating a velocity ratio may be helpful; e.g., the highest internal carotid artery PSV is divided by the PSV of the more distal common carotid artery.

**Stenosis of the Aortic Valve**

- A delay in the systolic upstroke is evident (Figure 2-24).
- Decreased peak systolic velocities (PSVs) are also usually seen. Therefore, PSVs may underestimate a stenosis.

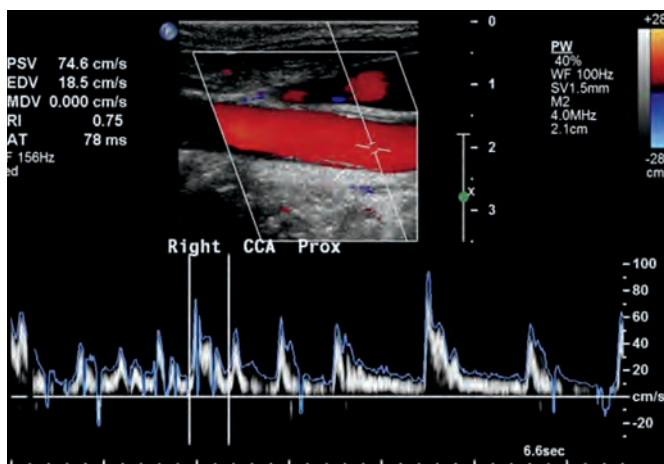


**Figure 2-23.** Severely irregular heart rhythm provides significant challenges in determining a reliable peak systolic velocity measurement.

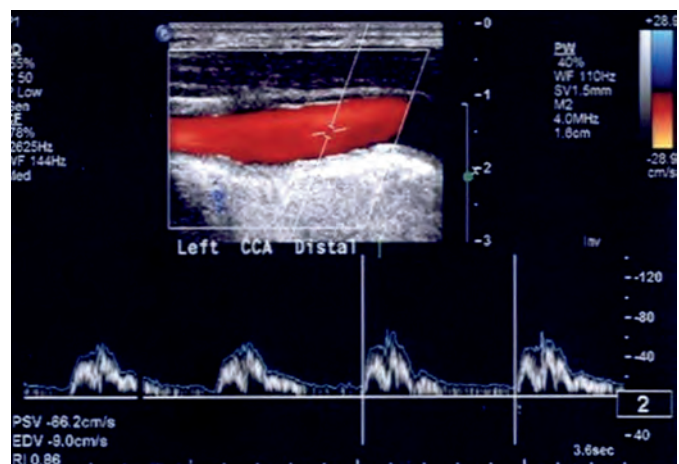


**Figure 2-24.** A delay in the systolic upstroke is evident with a stenosis of the aortic valve. A decrease in peak systolic velocities is also usually seen.


**Figure 2-25.** Double systolic peak (pulsus bisferiens) in the common carotid artery of a patient with aortic regurgitation/insufficiency.

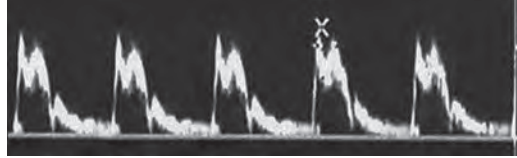

**Figure 2-26.** Pandsystolic flow reversal in the subclavian artery of a patient with aortic regurgitation/insufficiency.

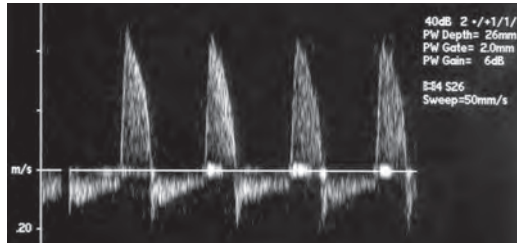

**Figure 2-27.** Significant flow reversal in the proximal aorta of a patient with severe aortic regurgitation/insufficiency.

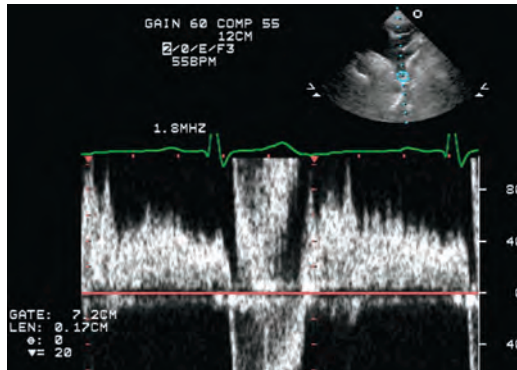

### Regurgitation/Insufficiency of the Aortic Valve


- A double systolic peak, sometimes referred to as *pulsus bisferiens*, may be evident (Figure 2-25). However, a similar finding, i.e., a second systolic peak, may be a normal finding in a young patient secondary to vessel wall compliance. In addition, one author\* suggests that there may be a relationship between a double systolic peak and aortic dissection.
- Diminished diastolic flow or even reversed flow throughout diastole may also be documented on spectral analysis (Figures 2-26 and 2-27).

### High Cardiac Output

- A systemic increase in peak systolic velocities may normally be evident in younger, athletic individuals or physically fit adults (Figure 2-28). A decrease in end-diastolic velocities may also be evident in a low-resistance vessel, e.g., the internal carotid artery. As previously mentioned, a double systolic peak may also be a normal finding in a young patient secondary to vessel wall compliance (Figure 2-25).

\*Burgess W: Recommended protocol for duplex ultrasound in common carotid artery dissection extending from the aortic arch. Paper presented at the Annual Conference of the Society for Vascular Ultrasound, San Francisco, California, May 2013.



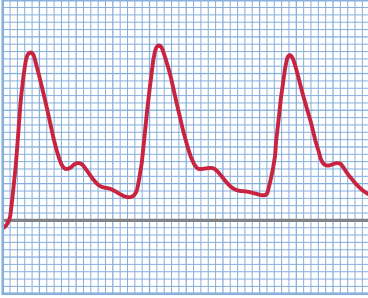

**Figure 2-40.** Because exercise produces a demand for blood to the muscles, a normally high-resistance arterial Doppler signal becomes low-resistance. In this analog waveform, the usual reversal below the baseline is seen as a forward reflection instead. This finding is quite normal after exercise.

## EFFECTS OF EXERCISE

Exercise should induce peripheral vasodilation in the microcirculation so that distal peripheral resistance diminishes and blood flow markedly increases.

Peripheral resistance—the resistance to blood flow caused by the ever-decreasing size of the vessels, especially in the microcirculation—changes in response to a variety of stimuli such as heat, cold, tobacco use, and emotional stress.

Vasoconstriction and vasodilation of the blood vessels within skeletal muscles are also influenced by sympathetic innervation fibers that function primarily to regulate body temperature.

Exercise is probably the best single vasodilator of high-resistance vessels within skeletal muscle.

Autoregulation also controls vasoconstriction and vasodilation. Autoregulation accounts for the ability of most vascular beds to maintain a constant level of blood flow over a wide range of perfusion pressures.


- Autoregulation does not function normally when perfusion pressure drops below a critical level.
- High-resistance vessels constrict in response to increased blood pressure and dilate in response to decreased blood pressure.

By decreasing resistance in the working muscle, exercise normally decreases reflection (flow reversal) of the Doppler flow signal in the exercising extremity.

**Example:** A low-resistance, monophasic Doppler flow signal (Figure 2-40) may be present normally in an extremity artery after vigorous exercise because the exercise causes peripheral dilatation and reduced flow resistance. This same low-resistance monophasic pattern is also seen pathologically when peripheral dilatation occurs in response to proximal arterial obstruction (see Figure 2-38).

On the other hand, a high-resistance signal (Figure 2-41) may occur from normal (physiologic) vasoconstriction at the arteriolar level or from distal arterial obstruction.

**Figure 2-41.** This high-resistance Doppler signal can occur with normal vasoconstriction at the arteriolar level. It can also occur proximal to distal arterial obstruction.

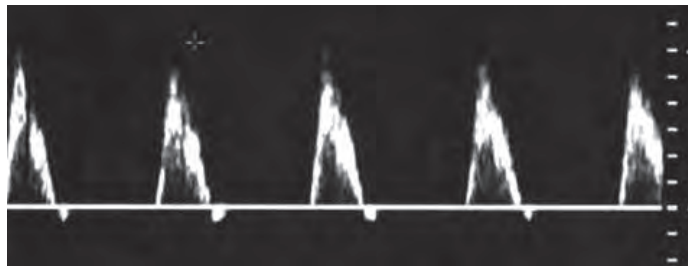

Additional notes:

- Proper characterization of velocity waveforms requires an understanding of both the normal flow characteristics of a particular artery and the physiologic status of the circulation supplied by the vessel.
- Questions need to be asked: Was the extremity cooled or warmed? Has it been exercised prior to the exam? Flow to a cool, vasoconstricted extremity will have pulsatile signals. Flow to a warm, vasodilated extremity will have continuous, steady signals.
- Proximal and distal pulsatility changes do not precisely differentiate between occlusion and severe stenosis.
- If good collateralization is present, proximal or distal Doppler velocity waveform qualities may not be altered.
- The distal effects of obstructive disease may be detectable in the presence of exercise or hyperemic evaluation.


## EFFECTS OF STENOSIS ON FLOW

Laminar flow review:


- Laminar flow has an even distribution of frequencies at systole, with the lower frequencies distributed at the walls (the boundary layer) and the higher frequencies in center stream.
- Stable flow through a relatively straight vessel is usually laminar, the layers of fluid slipping over one another with minimal friction normally.

A hemodynamically significant stenosis causes a major reduction in volume flow and pressure. A stenosis usually becomes hemodynamically significant when the cross-sectional area of the arterial lumen is reduced 75%, which corresponds to a diameter reduction of 50%.

- Diameter reduction is a one-dimensional measurement (Figure 2-42).
- Area reduction is a two-dimensional measurement (Figure 2-43).
- In hemodynamically significant stenoses, both pressure and flow volume decrease.

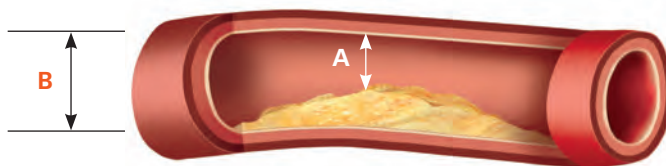

**Figure 2-42.** The nondiseased arterial segment is measured from wall to wall in a longitudinal approach (**B**) and then compared to the residual flow channel at the area of stenosis (**A**). A percentage diameter reduction is calculated.

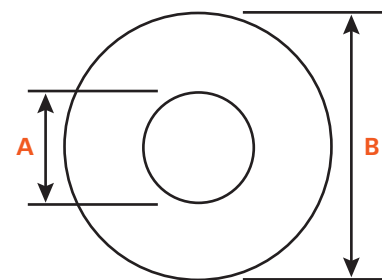

**Figure 2-43.** The original lumen (**B**) is measured and compared to the residual lumen (**A**) when the vessel is in a transverse approach. A percentage area reduction is calculated.

# Duplex Scanning and Color Flow Imaging of the Upper Extremities

*Capabilities and Limitations*

*Physical Principles*

*Technique*

*Interpretation*

*Brachial Artery Reactivity Testing*

## **CAPABILITIES AND LIMITATIONS**

### ***Capabilities***

- Localizes arterial stenosis or occlusion and evaluates degree of stenosis.
- Determines the presence or absence of aneurysm.
- Evaluates hemodialysis access graft or arterial bypass graft postoperatively.
- Detects arteriovenous fistulas or other unusual abnormalities.
- Evaluates medical treatment or surgery on a follow-up basis.

### Limitations

- Cannot be used in the presence of dressings, skin staples, sutures, or open wounds.
- Imaging is difficult around IV site.
- Hemodialysis access grafts:
  - ◆ Difficult to assess the anastomotic sites because of graft angulation.
  - ◆ Difficult to adequately evaluate the outflow vein secondary to increased collateral development.
- For diagnosis of Raynaud's phenomenon, previously described in Chapter 10 (page 113), physiologic studies combined with the clinical presentation (Chapter 3, pages 52–53) constitute the most useful diagnostic approach.

### PHYSICAL PRINCIPLES


- Duplex scanning combines real-time B-mode imaging (gray-scale evaluation) and Doppler spectral analysis (which is an analysis and display of the Doppler-shifted frequencies).
- Color flow imaging continues to provide the duplex information described above and in addition evaluates the Doppler flow information for its *phase* (direction toward or away from the transducer, on which basis color is assigned) and its *frequency content* (which determines the hue or shade of the assigned color).
- More in-depth coverage of this topic should be gleaned from textbooks focusing on physical principles of ultrasound scanning.

### TECHNIQUE

- The patient is supine with a small pillow under the head.
- The extremity is positioned close to the examiner.
- The arm to be evaluated is externally rotated and positioned at approximately a 45-degree angle from the body in what has been called the "pledge position."

For the patient with a dialysis graft, auscultate for a bruit and/or palpate for a "thrill" (vibration), bearing in mind that both graft stenosis and the high volume of blood flow through a patent dialysis graft will produce this effect.

- A multifrequency (e.g., 5–7 MHz) linear array transducer is used.
- The neck vessels are identified, with attention given to the innominate artery on the right. The left common carotid artery arises from the aortic arch.
- Duplex scanning (with or without color flow imaging) is performed at the following anatomic sites in this order:
  1. Subclavian artery
  2. Axillary artery

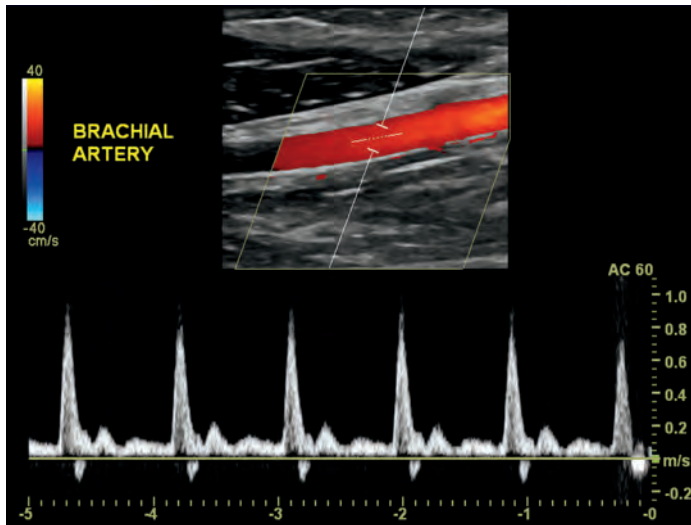

**Figure 12-1.** Spectral Doppler of the brachial artery, accompanied by color flow. Courtesy of Steve Bernhardt, BS, RDMS, RVT, RPhS.


3. Brachial artery
4. Radial artery
5. Ulnar artery
6. Palmar arch if necessary

Gray-scale imaging is used to observe vessel walls and to identify the presence of plaque and its morphology. Spectral analysis of these vessels, including evaluation of waveforms and measurement of peak systolic velocities, is used to assess blood flow. (See Figure 12-1.)


Color flow imaging is useful to discriminate flow deviations, flow channel narrowing within the vessel, and absence of flow where it should be, as well as assessing whether there is flow where there should not be (e.g., pseudo- or true aneurysm or inflamed lymph node).

- It is rather uncommon for arteries in the upper extremities to become stenotic. The main use for duplex and color flow imaging in the upper extremity is to evaluate hemodialysis access grafts, although other applications exist as well, including vein mapping and evaluation for thrombosis.
- Hemodialysis grafts are evaluated in the following order:
  1. Identify and evaluate inflow artery.
  2. Identify and evaluate arterial anastomosis.
  3. Identify and evaluate body of the graft.
  4. Observe for abnormalities, e.g., aneurysm, puncture site leaks, perigraft fluid collections.
  5. If color flow imaging is available, observe the image for frequency increases, turbulence, flow channel changes, and other deviations from normal.
  6. Identify and evaluate venous anastomoses.
  7. Identify and evaluate outflow vein.

**Figure 12-2.** Brescia-Cimino fistula being surgically created.

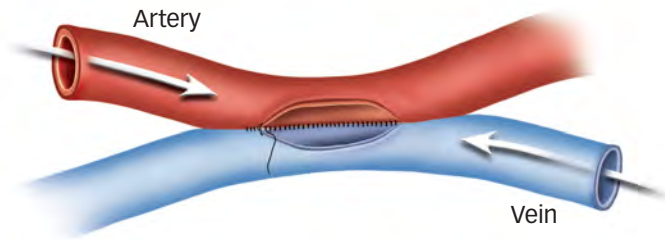

**Figure 12-3.** Straight arteriovenous dialysis access graft (brachial artery to axillary vein).

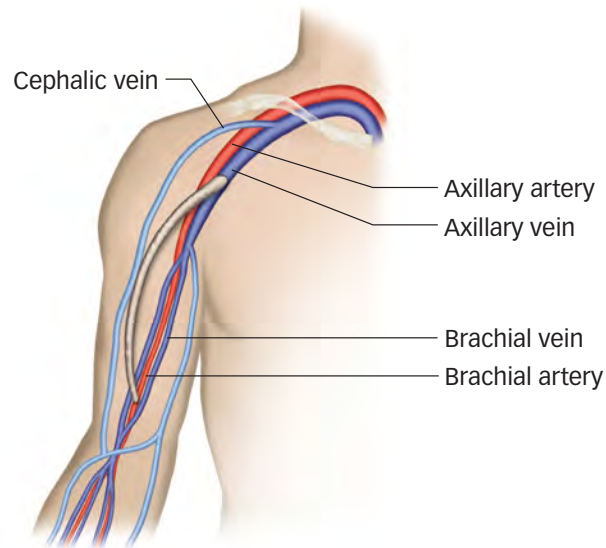


- Types of dialysis access include Brescia-Cimino fistulas (i.e., radial artery and cephalic vein) (Figure 12-2), straight synthetic grafts (see Figure 12-3), and looped synthetic grafts.
- Longitudinal and transverse approaches are used to evaluate the gray-scale image for thrombus, wall thickening, and other changes.
- Doppler peak systolic velocities (PSVs) are measured at the aforementioned sites (depending on which study is being performed) and also at other locations of interest (e.g., for stenosis, preocclusion).
- Documentation consists of storing the digital images (e.g., via the Picture Archiving and Communication System [PACS]) or making hard-copy prints.

## **INTERPRETATION**

### ***Stenosis***

- Currently, there are no criteria for classifying upper extremity disease as there are for the lower extremities.
- Peak systolic velocities vary widely with changes in skin temperature. The Doppler signal quality is usually triphasic (Figure 12-4A), but when the hand is quite warm or the arm has been exercised, the Doppler signal quality becomes lower-resistance (continuous through diastole). (See Figure 12-4B.) The effects of cooling were previously discussed in Chapter 10, page 111.

**Figure 14-2.** In this transverse image of an abdominal aortic aneurysm, at least 7 cm in diameter, the A/P measured distance is depicted with black X's, while the lateral (also referred to as transverse) measured distance is depicted with white X's. It contains a mixture of echoes that are more heterogeneous than homogeneous. Peak systolic velocities vary according to the extent of thrombus. Large lumen (less thrombus) = lower velocities; small lumen (more thrombus) = higher velocities; preocclusive = lower velocities.

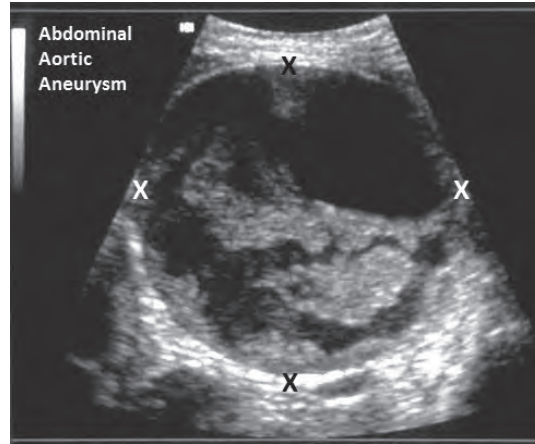


- ◆ The examiner should note the type of aneurysm—true (i.e., fusiform or saccular), false, or dissecting—and the presence of thrombus, if any. See the example of a thrombotic aneurysm in Figure 14-2.
- ◆ The most frequent complication and danger of an AAA is rupture, but it can also embolize. Both thrombosis and embolization are considered primary complications of peripheral arterial aneurysms. It is not uncommon for both abdominal and peripheral aneurysms to contain varying amounts of thrombus.


### Renal Artery and Kidney—Technique

See “General Remarks” on page 140. It is helpful to know that many of the patients undergoing this study present with hypertension (controlled or not well controlled) and that many of these hypertensive patients have *renovascular* hypertension. Renovascular hypertension, a secondary form of high blood pressure, is often caused by renal artery stenosis (which can be secondary to atherosclerosis or fibromuscular dysplasia) or occlusion. The narrowing of the artery reduces blood flow to the kidney and, in response, the kidney produces the enzyme renin. Through a series of complicated interactions, there is a conversion of angiotensinogen to angiotensin II, which results in constriction of the blood vessels. Through another process, sodium and fluid retention also occur.


- The examiner begins in the longitudinal plane, acquiring velocity data from the celiac and superior mesenteric arteries.
- The peak systolic velocity of the aorta is obtained just distal to the superior mesenteric artery.
- Using a transverse approach, the examiner locates the renal arteries. The left renal vein is a good landmark for identifying the left renal artery, as depicted in Figure 14-3. See also Figures 14-4, 14-5, and 14-6, which show not only the left renal vein but also the renal artery branches.
- Kidney size and morphology are evaluated bilaterally. Figure 14-7 shows a normal-sized kidney with color flow clearly demonstrating the segmental arteries, as well as the renal artery.

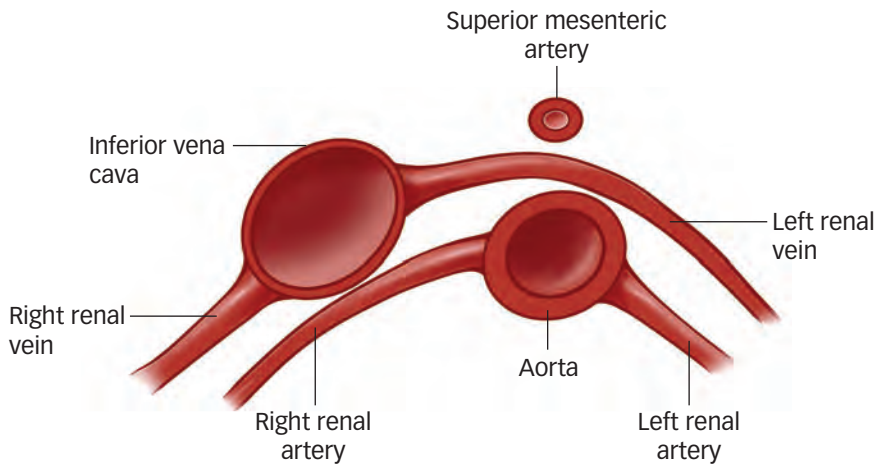


**Figure 14-3.** Transverse view of the abdomen showing the position of the left renal vein in relation to the left renal artery.

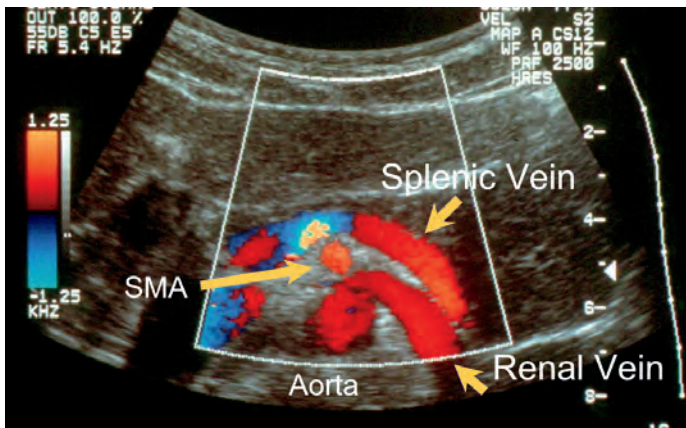


**Figure 14-4.** Transverse view of the abdominal aorta also showing the superior mesenteric artery (SMA), left renal vein, and splenic vein. This image is a good reminder to watch the color bars for flow direction. The veins, in this example, are not blue.

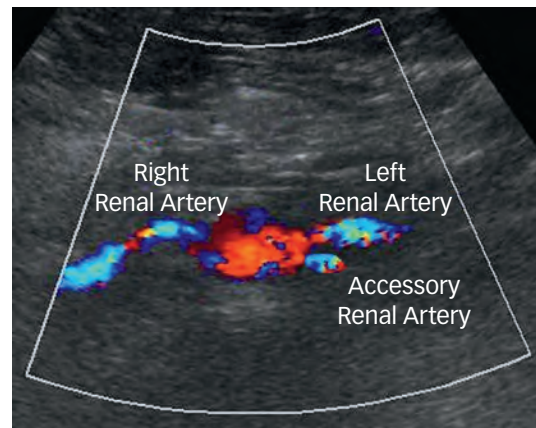

**Figure 14-5.** Transverse image of the abdominal aorta, showing two left renal arteries and a single right renal artery.

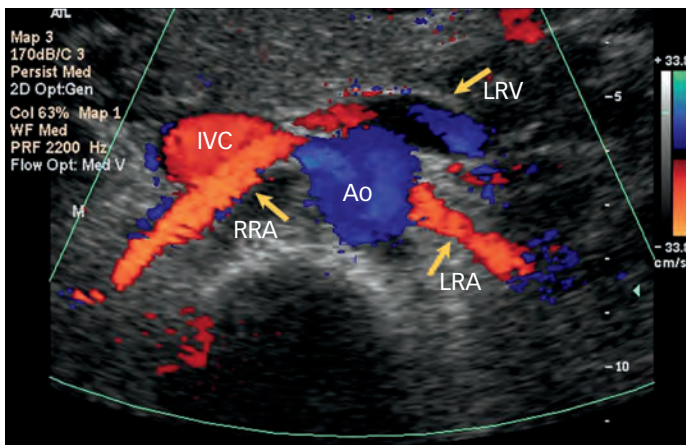


**Figure 14-6.** Transverse image of the abdominal aorta (Ao), showing the following vessels with an accompanying arrow: right (RRA) and left (LRA) renal arteries and left renal vein (LRV), which can be insonated with an appropriate Doppler angle. Note also the inferior vena cava (IVC), to the left of the aorta and the left renal vein. Courtesy of UC Davis Medical Center.

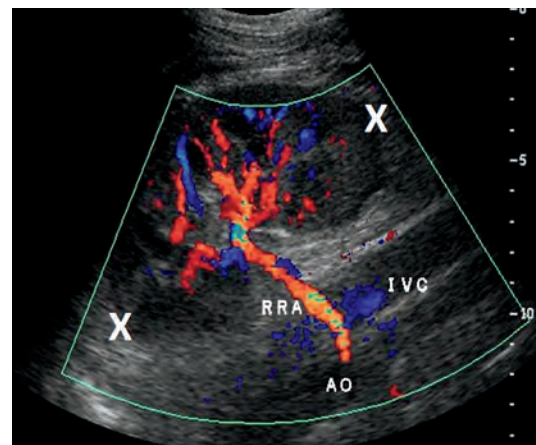


**Figure 14-7.** Normal-sized kidney with color flow clearly showing the segmental arteries, as well as the right renal artery (RRA) depicted from its origin at the aorta (Ao) into the kidney. The X's shown on the ends of the kidney depict the start and end points for length measurement. Courtesy of UC Davis Medical Center.



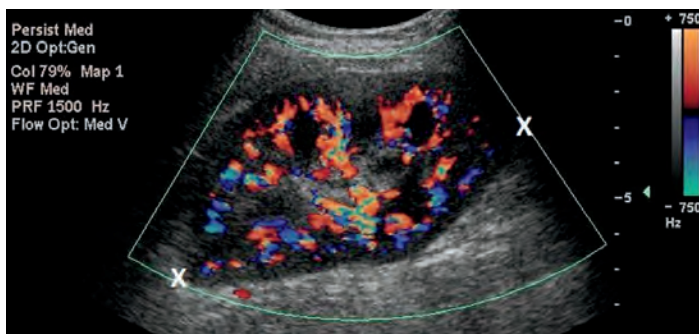

**Figure 14-13.** Longitudinal view of a normal kidney. The X's located at the ends of the kidney show how the kidney length is measured.

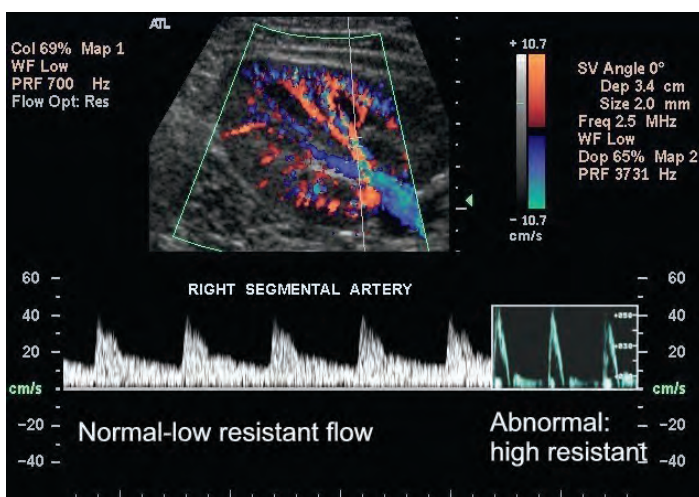

**Figure 14-14.** Transverse view of the right kidney, using a coronal approach. Doppler signals are obtained from a segmental artery within the kidney. Normally, the signals are low-resistance, as noted on the image with the first set of waveforms. In an abnormal (e.g., nephrosclerotic) kidney, the Doppler signal is high-resistance, as seen to the right.

The results are interpreted as follows:

- ◆ *Normal:*  $\geq 0.2$
- ◆ *Abnormal:*  $< 0.2$

An EDR or PRR of  $< 0.2$  indicates an increase in resistance within the kidney parenchyma. An example of the type of waveform obtained in this situation appears in Figure 14-14.\*

- Another ratio used (to a lesser degree) to determine whether resistance is increasing in the kidney (especially the transplanted kidney) is Pourcelot's ratio/resistivity index (RI), which is calculated as follows:

$$\frac{PSV - EDV}{PSV}$$

- ◆ *Normal:*  $< 0.7$
- ◆ *Abnormal (increased resistance):*  $\geq 0.7^{**}$


\*Neumyer MM: Native renal artery and kidney parenchymal blood flow duplex evaluation with/without color flow imaging. In Rumwell CB, McPharlin MM, Strandness DE, et al (eds): *Vascular Laboratory Policies & Procedures Manual*. Pasadena, Davies Publishing, 2008.

\*\*Rifkin MD, Needleman L, Pasto ME, et al: Evaluation of renal transplant rejection by duplex Doppler examination: value of the resistive index. *Am J Roentgenol* 148:759–762, 1987.

Example of velocity findings in a kidney with abnormal arterial flow patterns:

- ◆ *Doppler findings*: PSV = 70, EDV = 10
  - ◆ *End-diastolic ratio (EDR)*:  $10/70 = 0.14$
  - ◆ *Resistivity index (RI)*:  $\frac{70 - 10}{70} = 0.86$
- The renal resistive-index value (RRIV)\* also uses the peak systolic velocity and the end-diastolic velocity obtained from the kidneys' segmental arteries. The purpose of using this value is to help prospectively identify patients whose renal function or blood pressure will improve after correction of the renal artery stenosis. A lower-resistance RRIV value of  $<0.75$  is associated with improvement in both blood pressure and renal function after the correction of renal artery stenosis. Conversely, an RRIV of  $\geq 0.75$  is a strong predictor of worsening renal function and lack of blood pressure improvement, despite the correction of renal artery stenosis. Calculation of the RRIV is as follows:  $1 - \text{EDV/PSV}$ . When you compare this to the EDR or PRR, described above, note that the sum of the PRR and the RRIV is 1; that is, subtract the PRR from 1 and you will get the RRIV, and vice versa.
  - Proximal high-grade stenosis or occlusion of the renal artery may result in a dampened (prolonged-upstroke) but still low-resistance waveform (see Figure 14-15A). The angiographic result (see Figure 14-15B) is an example of what can cause the poor flow depicted in Figure 14-15A.
  - Some references describe the usefulness of evaluating the kidney arterial flow by determining the acceleration time (AT) and acceleration index (AI).
    - ◆ The AT is the time interval from the onset of systole to the initial peak and is reported in milliseconds (msec). A proximal stenosis of  $\geq 60\%$  diameter reduction is most likely to produce an AT of  $\geq 100$  msec.
    - ◆ The AI describes the slope of the Doppler velocity waveform. It is calculated as the change in velocity between the onset of systole and the

**Figure 14-15. A** Spectral Doppler arterial signal obtained from the kidney in the presence of a proximal high-grade stenosis or occlusion of the renal artery. This waveform pattern is termed *tardus parvus*. **B** Angiogram showing an occlusion (arrow) of the left renal artery. Collateral vessels are providing poor-quality arterial flow to the kidney.

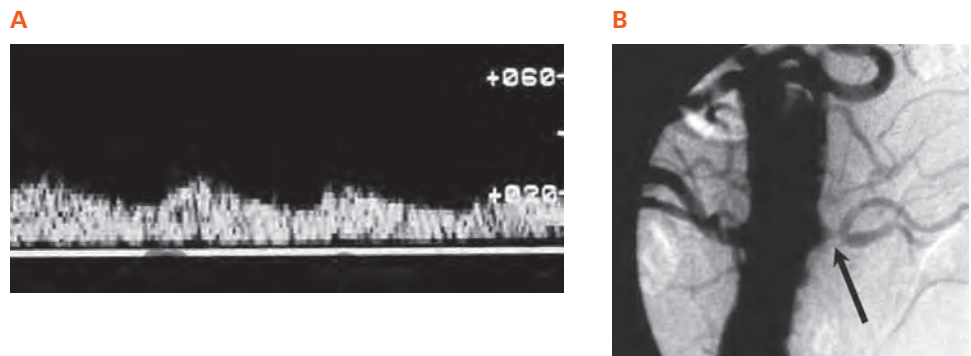


\*Yuksel UC, Anabtawi AG, Cam A, et al: Predictive value of renal resistive index in percutaneous renal interventions for atherosclerotic renal artery stenosis. *J Invasive Cardiol* 24:504–509, 2012.


systolic peak (cm/sec) divided by the AT. The units of measure for AI are  $\text{cm}/\text{sec}^2$ . A positive result (consistent with a  $\geq 60\%$  diameter reduction) is defined as  $\text{AI} \leq 291 \text{ cm}/\text{sec}^2$ .\*, \*\*

### Mesenteric Arteries—Technique

- Patient history: Patients who present with a history of dull, achy, or crampy abdominal pain 15–30 minutes after meals may suffer from *mesenteric ischemia*. Mesenteric ischemic pain is also known as *mesenteric angina*. This condition may be caused by a stenosis or occlusion of vessels such as the superior mesenteric, celiac, or inferior mesenteric arteries.
- Mesenteric ischemia may be acute or chronic. Both conditions are difficult to diagnose. The study described here assists with the diagnostic process, but an arteriogram is essential for diagnosis.
- It is essential that this study always be performed on a fasting patient.
- Peak systolic and end-diastolic velocities are obtained in the longitudinal plane from the following arteries:
  - ◆ Celiac trunk (artery), hepatic artery, and splenic artery, all of which often require a transverse approach to the aorta (see Figure 14-16)
  - ◆ Proximal, mid, and distal superior mesenteric artery (see Figure 14-17)

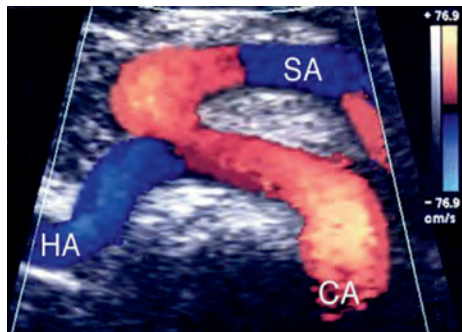

**Figure 14-16.** Celiac artery (CA)—i.e., celiac trunk—with branches. The hepatic artery (HA) takes blood to the liver; the splenic artery (SA) takes blood to the spleen. All Doppler flow is low-resistance because organs have a constant need for blood supply. Courtesy of Steve Bernhardt, BS, RDMS, RVT, RPhS.

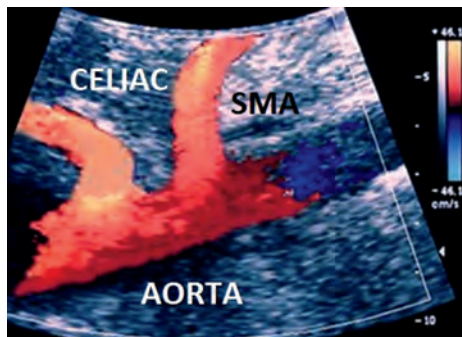

**Figure 14-17.** Color flow Doppler image using a low-frequency transducer in the longitudinal (sagittal) approach demonstrates the celiac artery followed by the superior mesenteric artery (SMA) branching off the aorta. Courtesy of Steve Bernhardt, BS, RDMS, RVT, RPhS.

\*Issacson JA, Zierler RE, Spittell PC, et al: Noninvasive screening for renal artery stenosis: comparison of renal artery and renal hilar duplex scanning. *J Vasc Technol* 19:105–110, 1995.


\*\*Martin RL, Nanra RS, Włodarczyk J, et al: Renal hilar Doppler analysis in the detection of renal artery stenosis. *J Vasc Technol* 15:173–180, 1991.

**Table 14-1** Summary criteria for interpreting pre- and postprandial findings in the superior mesenteric and celiac arteries.


Description	SMA	Celiac Artery
<i>Preprandial (fasting):</i>		
PSV	High	High
EDV	Low	High
Flow reversal	Yes	No
<i>Postprandial (food challenge):</i>		
PSV	Marked increase	No change
EDV	Marked increase	No change
Loss of flow reversal	Yes	N/A
<i>Velocity criteria:</i>		
Normal PSV*	110–177 cm/sec	50–160 cm/sec
<i>Stenosis criteria*</i>		
(includes poststenotic turbulence)	PSV $\geq$ 275 cm/sec predicts $\geq$ 70% diameter reduction	PSV $\geq$ 200 cm/sec predicts $\geq$ 70% diameter reduction

\*Criteria and other information from various publications by Gregory L. Moneta, MD, Oregon Health & Science University Hospital, Portland, Oregon.

- Celiac band syndrome, which is an extrinsic compression of the celiac artery origin by the median arcuate ligament of the diaphragm, is a fairly frequent anatomic finding, especially among young, athletic women. A reversible celiac artery stenosis occurs and is rarely the cause of clinical symptoms. In such cases the high-velocity signals indicative of stenosis are substantially improved with deep inspiration and return with expiration. See Figure 14-23. An audible bruit can be auscultated, and during duplex/color flow imaging a “color bruit” may appear during exhalation.
- The celiac artery branches into the hepatic and splenic arteries (see Figure 14-16). The liver and spleen have fixed metabolic requirements and are not likely to be influenced by the postprandial state. Therefore, high peak systolic and end-diastolic velocities are normally evident in these vessels at all times.

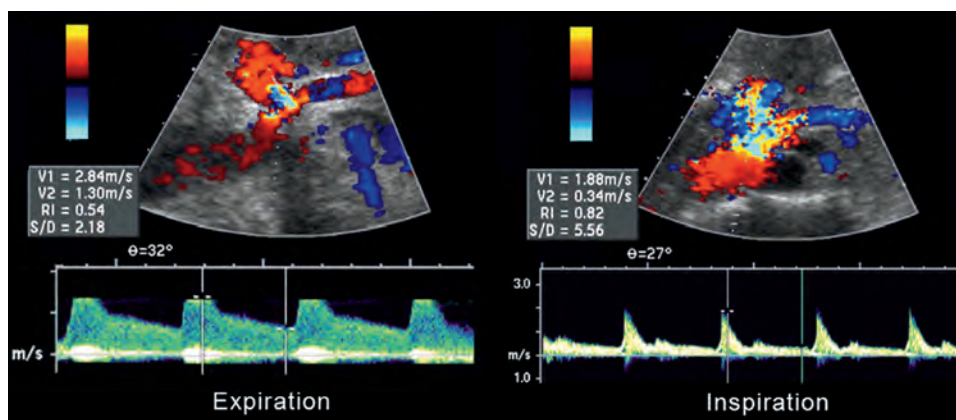

**Figure 14-23.** Longitudinal image of the abdominal aorta and celiac artery. The Doppler spectral waveforms were recorded at both phases of respiration. Extrinsic compression of the celiac artery during expiration is a frequent finding. In this example, it causes peak systolic velocities and end-diastolic velocities to increase ( $>284/130$  cm/sec) with aliasing noted. During inspiration, velocities decrease to 188/34 cm/sec. Courtesy of Steve Bernhardt, BS, RDMS, RVT, RPhS.

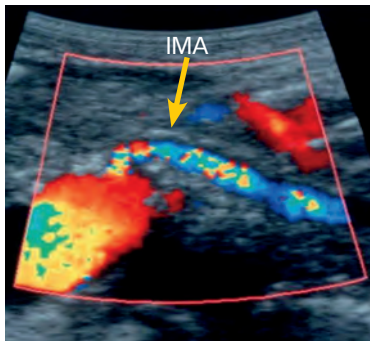


**Figure 14-24.** Inferior mesenteric artery (IMA) easily observed in this image because of a superior mesenteric artery (SMA) occlusion (not seen on this image).

- The hepatic artery and splenic artery Doppler signals should be low-resistance in quality because each provides arterial inflow to a particular organ. There are no specific peak systolic velocity criteria for what suggests stenosis in either of these vessels (other than the stenosis profile described on pages 40–41), but if the celiac artery is occluded, retrograde flow may occur in the hepatic artery and, if so, should be identified during the duplex evaluation.
- Because of its small caliber, it is quite difficult to locate the inferior mesenteric artery (IMA), which branches off the distal aorta. Easy detection and ultrasonographic dominance of the inferior mesenteric artery (see Figure 14-24) may suggest occlusion of the superior mesenteric artery.
- Of interest is the ability of the IMA to provide collateralization in response to an occlusion of the SMA. Two possible collateral connections between the SMA and IMA include the arc of Riolan and the marginal artery of the colon (also known as the marginal artery of Drummond).
- Chronic mesenteric ischemia is suggested when there are abnormal findings for at least two of the three mesenteric vessels (celiac, superior mesenteric, and inferior mesenteric arteries) in symptomatic patients.


## **ORGAN TRANSPLANTS**

### ***Liver Transplant (Allograft)***

- Liver transplantation is becoming a more frequently utilized procedure to treat patients with end-stage liver disease. In the pre- and postoperative assessment of these patients, the examiner uses a multifrequency (e.g., 3–5 MHz) linear or phased array transducer.
- The preoperative duplex evaluation of candidates for liver transplantation includes documenting the patency of the portal vein, splenic vein, superior mesenteric vein, hepatic veins, inferior vena cava, and hepatic artery while observing for abnormalities (i.e., tumors in the liver and surrounding tissue) and/or other vasculature (see Figure 14-25). Doppler waveforms are evaluated according to the appropriate criteria. It is also important to determine the status of the biliary tree. Abnormalities such as thrombosis of the portal vein can make transplantation difficult, if not impossible.
- Postoperatively, duplex ultrasonography is used to document patency of the portal vein (Figure 14-26), splenic vein, superior mesenteric vein, hepatic veins, inferior vena cava, and hepatic artery, as well as portal vein flow direction and vessel size. Normal portal vein size should be  $\leq 1$  up to 1.5 cm.
- Postoperative complications include allograft rejection, pseudoaneurysm, hepatic infarction, and thrombosis of the portal vein, inferior vena cava, and/or hepatic artery. Obviously, the hepatic artery inflow becomes crucial to the viability of the liver.
  - ◆ The location of this artery provides a challenge to obtaining a proper Doppler angle. If the Doppler angle is high (i.e., 80–90 degrees), no flow

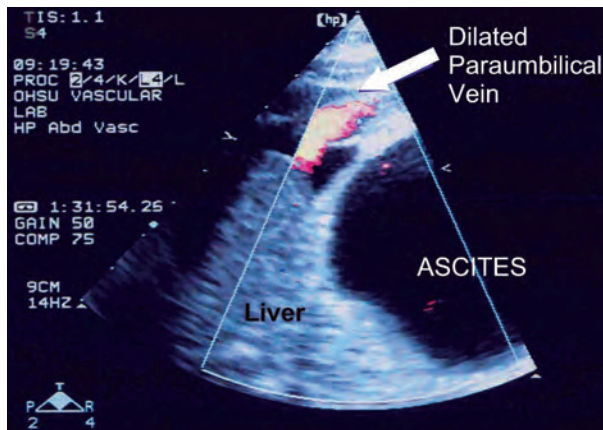

**Figure 14-25.** Sagittal/oblique view of the liver revealing ascites and a dilated paraumbilical vein. It is not unusual to observe one or both of these findings in patients being evaluated for liver transplantation. A sudden increase of ascites often precedes portal vein thrombosis, especially in cirrhotic patients. Dilatation of the paraumbilical vein is a specific finding for portal hypertension, in which normally hepatopetal flow becomes hepatofugal as blood flow is diverted to the systemic venous circulation through a collateral pathway.

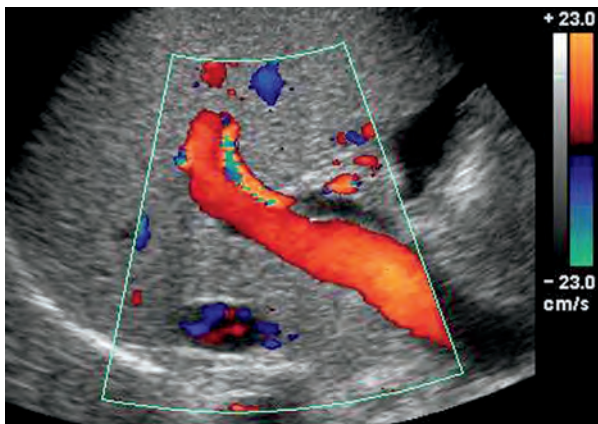

**Figure 14-26.** Duplex scan of the liver with a patent, normal portal vein. Flow direction (red color) is toward the liver (hepatopetal).

may be detected because of the poor or nonexistent Doppler shift. Before concluding that the hepatic artery is occluded, the examiner must be certain that an appropriate Doppler angle has been used. If possible, the proximal hepatic artery is examined where it branches off the celiac artery to determine the presence or absence of flow.

- ◆ Acute rejection will cause liver dysfunction. It has been thought by some that rejection would cause arterial resistance (measured by the resistivity index) to increase. Others, however, suggest that rejection is a cellular process that does not affect the peripheral vessels.

### **Renal Transplant (Allograft)**

- Both living related–donor and cadaveric kidneys are implanted most commonly in the right iliac fossa.
- Postoperative examination is facilitated if the examiner knows what type of anastomosis was used. The transplanted renal artery vessels are anastomosed either end-to-side with the external iliac artery or end-to-end with the internal iliac artery. The venous anastomosis connects the renal vein to the external iliac vein in an end-to-side fashion. Most often, the donor ureter is implanted into the bladder directly.

# Vascular Technology

AN ILLUSTRATED REVIEW

5th EDITION

Claudia Rumwell, RN, RVT, FSVU | Michalene McPharlin, RN, RVT, RVS, FSVU

For registry candidates, technologists in training, cross-training sonographers, clinical personnel, and RPVI candidates, here is the completely revised, updated 5th edition of the best and most trusted vascular technology review text available—now, for the first time, in full color. Written and peer-reviewed by nationally renowned authors and vascular experts, it remains a silver-bullet review of everything a registry candidate needs to know to pass the ARDMS specialty exam in vascular technology. Also used in labs and diagnostic ultrasound programs across the country, it is perhaps the most widely read vascular reference in print. This new edition includes more than 500 color illustrations and images, photo-documented physiologic exam techniques, 120 self-assessment case studies, and an integral application for 15 SDMS-approved CME credits. RVT, RPVI, and RVS candidates and DMS students can combine this Step 1 review text with Davies' Step 2 mock exam *Vascular Technology Review* and Step 3 *ScoreCards for Vascular Technology* for powerful 1-2-3 Step Ultrasound Education & Test Preparation. Ready to score? With this proven review, you can.

## Other publications of interest

***Vascular Technology Review***, edited by Donald P. Ridgway, RVT, and D.E. Strandness, Jr., MD. This completely revised edition illuminates the facts you need to know to pass the Vascular Technology specialty exam, hones your test-taking skills, and reveals your strengths and weaknesses by exam topic. Based on the ARDMS exam outline, it delivers nearly 600 Q&A items with explanations, references, and lots of images. 12 hours CME credit.



***Vascular Physics Review***, edited by Barton A. Bean, RVT. Here's the continuously updated vascular physics review RVT candidates rely on. Approximately 500 illustrated question/answer/explanation items in registry format to simulate the vascular-specific physics questions on the vascular technology and PVI examinations. 7.5 hours CME credit.

***ScoreCards for Vascular Technology***, by Cindy Owen, RDMS, RVT, and D.E. Strandness, Jr., MD. This portable flip-card study system exercises your ability to think fast, recall key facts, and apply knowledge, wherever you are. And it's fun. More than 400 Q&A items keyed to the new registry outline, 50 image-based questions in the Image Gallery, explanations, references. 7.5 hours CME credit.



***Vascular Anatomy and Physiology***, by Ann C. Belanger, RN, RVT. The perfect introduction to vascular anatomy and physiology for everyone who seeks a clear and simple presentation of the facts they must know, including hemodynamics. Highly recommended. 5 hours CME credit.

***Introduction to Vascular Scanning***, 4th edition, by Donald P. Ridgway, RVT. For novice scanners and for sonographers and echocardiographers cross-training in vascular ultrasound, here is the new and improved version of Don Ridgway's very popular, unabashedly practical, and famously unique guide to performing vascular studies—now with new chapters on the Doppler principle, Those Darn Doppler Angles, and other vascular diagnostic modalities. 18 hours CME credit.



## CD-ROM Mock Exams.

Powerful and fun, these interactive wonders simulate the ARDMS exams

and deliver content you can trust in Learn Mode or Test Mode. Each includes hundreds of questions in registry format with answers and clear explanations, along with scores of image-based questions, references, tutorials, and automatic performance analysis. A special bonus: 15 hours each of CME credit for Physics, Abdomen, Ob/Gyn, and Vascular; 12 for Fetal Echo; and 10 for Breast.



Claudia Rumwell, RN, RVT, FSVU

Claudia Rumwell lectures on registry preparation, provides consulting services, and continues to write. She is a past member of the board of the Society for Vascular Ultrasound, is a recipient of the SVU Distinguished Service Award, and has been honored with the status of SVU Fellow. She previously served in the following roles for the Division of Vascular Surgery at Oregon Health & Science University: Technical Director of the Vascular Lab, Clinic Nurse Coordinator, Research Nurse Coordinator, and Instructor in Surgery.



Michalene McPharlin, RN, RVT, RVS, FSVU

Micky McPharlin lectures widely on registry preparation and vascular topics, as well as providing consulting services. She also serves on the editorial board of the *Journal for Vascular Ultrasound*. Like Claudia, she is a recipient of the Society for Vascular Ultrasound Distinguished Service Award and has been honored with the status of SVU Fellow. She is the former Technical Director of the Vascular Laboratory of Henry Ford Hospital in Detroit, Michigan, and the former Program Director of the Vascular Ultrasound Program at Baker College of Auburn Hills, Michigan.

ISBN 978-0-941022-85-9




www.daviespublishing.com



# DAVIES

*Introduction to*

# Vascular Ultrasonography

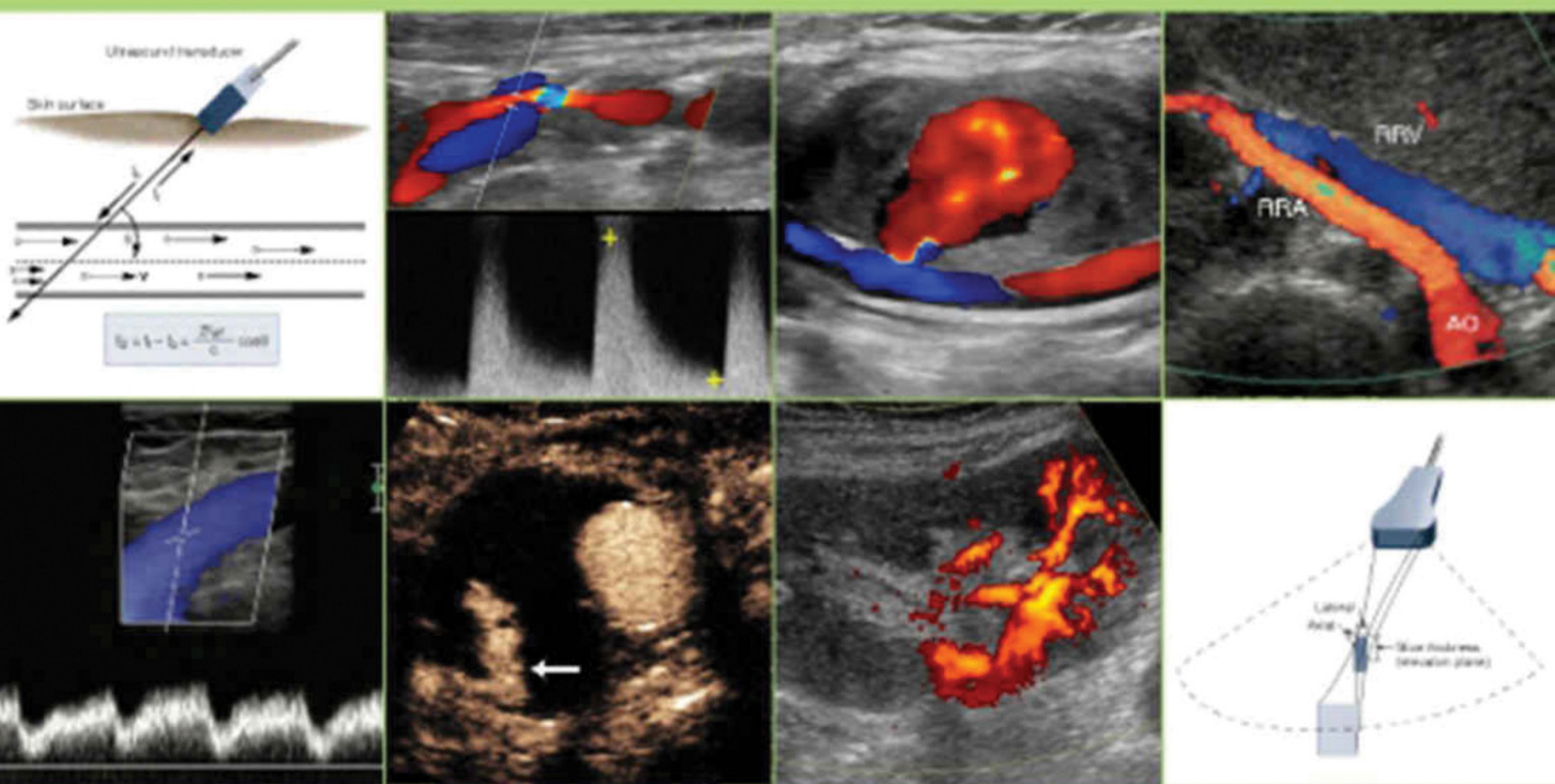

Pellerito  
Polak



# Introduction to Vascular Ultrasonography

---


7th Edition

**John S. Pellerito, MD, FACR, FSRU, FAIUM**

Professor of Radiology  
Department of Radiology  
Zucker School of Medicine at Hofstra/Northwell  
Hempstead, New York  
Vice Chairman  
Department of Radiology  
Northwell Health System  
Manhasset, New York

**Joseph F. Polak, MD, MPH, FACR, FAIUM**

Professor of Radiology  
Department of Radiology  
Tufts University School of Medicine  
Director  
Ultrasound Reading Center  
Boston, Massachusetts  
Chief of Radiology  
Department of Radiology  
Lemuel Shattuck Hospital  
Jamaica Plain, Massachusetts

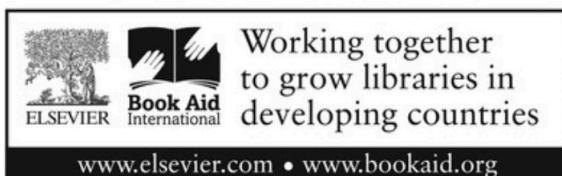

# Copyright

---

Elsevier  
1600 John F. Kennedy Blvd.  
Ste 1800  
Philadelphia, PA 19103-2899

INTRODUCTION TO VASCULAR ULTRASONOGRAPHY, SEVENTH EDITION ISBN: 978-0-323-42882-8  
**Copyright © 2020 by Elsevier, Inc. All rights reserved.**

No part of this publication may be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopying, recording, or any information storage and retrieval system, without permission in writing from the publisher. Details on how to seek permission, further information about the Publisher's permissions policies and our arrangements with organizations such as the Copyright Clearance Center and the Copyright Licensing Agency, can be found at our website: [www.elsevier.com/permissions](http://www.elsevier.com/permissions).

This book and the individual contributions contained in it are protected under copyright by the Publisher (other than as may be noted herein).

## Notice

Practitioners and researchers must always rely on their own experience and knowledge in evaluating and using any information, methods, compounds or experiments described herein. Because of rapid advances in the medical sciences, in particular, independent verification of diagnoses and drug dosages should be made. To the fullest extent of the law, no responsibility is assumed by Elsevier, authors, editors or contributors for any injury and/or damage to persons or property as a matter of products liability, negligence or otherwise, or from any use or operation of any methods, products, instructions, or ideas contained in the material herein.


Previous editions copyrighted 2012, 2005, 2000, 1992, 1986, and 1983.

**Library of Congress Control Number: 2019949237**

*Content Strategist:* Russell Gabbedy  
*Senior Content Development Specialist:* Joanie Milnes  
*Publishing Services Manager:* Catherine Jackson  
*Senior Project Manager:* Claire Kramer  
*Design Direction:* Amy Buxton

Printed in Canada

Last digit is the print number: 9 8 7 6 5 4 3 2 1

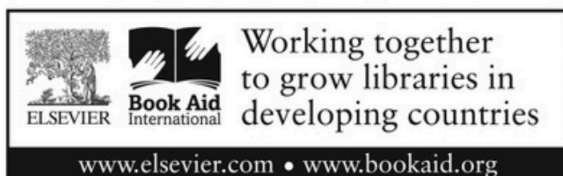

# Preface

---

The seventh edition of *Introduction to Vascular Ultrasonography* is a significant update to our previous editions. My co-editor, Joseph F. Polak, and I have pulled together the topics and authors we thought represented the best of vascular ultrasound imaging. Having collaborated together for more than two decades with our vascular course and the sixth edition of this book, we feel that this edition has met our goal to present the definitive text in vascular ultrasound. Jo and I have contributed to 22 of the 35 chapters. Not only have all the chapters in this edition been revised to deliver the newest techniques, protocols, and topics in vascular ultrasound, we also invited several new experts to provide their perspectives and experience in significant areas. For example, the physics section has been completely revamped by Dr. Fred Kremkau, the favored speaker and authority in ultrasound physics. We welcome Heather Gornik, MD, the current president of IAC vascular testing, for her contribution on credentialing, accreditation, and quality in the vascular laboratory. We also are excited to have several world-renowned experts provide us the state of the art in ultrasound contrast applications for vascular imaging. Of course, we brought back many of our favorite authors who received positive reviews from our previous edition.

In addition to multiple new figures and illustrations, a major enhancement to this version of *Introduction to Vascular Ultrasonography* is the integration of Practical Tips sections throughout all the chapters. These Tips are intended to focus the reader on major teaching points and pearls for the successful performance and interpretation of vascular studies.

We are extremely proud to present the seventh edition of *Introduction to Vascular Ultrasonography*. We hope that this book will provide guidance to students, technologists, sonographers, and all practitioners of vascular ultrasound to improve patient diagnosis and management.

*John S. Pellerito, MD, FACR, FSRU, FAIUM*

# Table of Contents

---

## Section 1—Basics

### 1 The Hemodynamics of Vascular Disease

- Abstract
- Overview
- Introduction
- Physiologic Factors Governing Blood Flow and Its Characteristics
- Effects of Arterial Obstruction
- Venous Hemodynamics
- Summary
- References

### 2 Principles and Instruments of Ultrasonography

- Introduction
- Ultrasound Principles
- Transducers
- Sonographic Instruments
- Advanced Features
- Doppler Principles
- References

### 3 Doppler Flow Imaging and Spectral Analysis

- Introduction
- Blood Flow Patterns
- Listening to the Auditory Frequency Spectrum
- What You Need to Know About Waveform Analysis
- Diagnosis of Arterial Obstruction
- Color Flow Ultrasound Imaging
- Three-Dimensional Vascular Imaging
- B-Mode Flow Imaging
- References

## Section 2—Cerebral Vessels

### 4 Anatomy of the Cerebral Structure

- Abstract
- Introduction
- Vascular Anatomy
- Cerebral Hemodynamics
- Collateralization
- Summary
- References

### 5 Carotid Sonography

- Abstract
- B-Mode Ultrasound of the Normal Carotid Artery Wall
- Normal Blood Flow Characteristics
- Vessel Identity
- Examination Protocol
- The Examination Sequence
- References

### 6 Evaluating Carotid Plaque and Carotid Intima-Media Thickness

- Abstract
- Introduction
- Pathogenesis of Carotid Atherosclerosis
- Imaging Protocol: IMT and Plaque
- Plaque Characterization
- Plaque Classification Schemes
- Plaque Surface Features
- Plaque Neovascularity
- Perspective on Plaque Characterization
- Paradigm Shift: Intima-Media Thickness
- Carotid Plaques: Extent, Severity, and Follow-Up
- References

### 7 Ultrasound Assessment of Carotid Stenosis

- Abstract
- Introduction
- Clinical Background
- Imaging of Stenosis
- Diagnostic Criteria
- Plaque Characterization
- Technical Considerations and Pitfalls
- The Common Carotid Artery
- The External Carotid Artery
- The Innominate Artery
- Summary
- References

### 8 How to Assess Difficult and Uncommon Carotid Cases

- Introduction
- Carotid Occlusion
- Carotid Artery Dissection

Carotid Pseudoaneurysm  
Carotid Arteriovenous Fistula  
Fibromuscular Dysplasia  
Carotid Body Tumor (Glomus Tumor)  
Difficult Carotid Cases

## 9 Ultrasound Assessment of the Vertebral Arteries

Abstract  
Introduction  
Examination Technique  
Vertebral Artery Hemodynamics: Qualitative Assessment  
Vertebral Artery Hemodynamics: Quantitative Assessment  
Correlative Imaging With Magnetic Resonance Angiography  
and Computed Tomographic Angiography  
Treatment of Disease  
Summary  
References

## 10 Ultrasound Assessment of the Intracranial Arteries

Introduction  
Examination Techniques  
Diagnostic Parameters for Specific Clinical Applications  
Summary  
References

## Section 3—Extremity Arteries

### 11 Anatomy of the Upper and Lower Extremity Arteries

Abstract  
Introduction  
Upper Extremity  
Lower Extremity  
References

### 12 Physiologic Testing of Lower Extremity Arterial Disease

Abstract  
Introduction  
Instrumentation  
Physiologic Testing Procedures  
Diagnostic Algorithm  
Clinical Applications  
Summary  
References

### 13 Assessment of Upper Extremity Arterial Disease

Abstract  
Introduction  
Basic Anatomy  
Obstructive Arterial Disease  
Thoracic Outlet Evaluation  
Digital Artery Evaluation  
Raynaud Disease/Phenomenon  
Color Doppler Imaging and Duplex Ultrasound of Vasculitis  
Ultrasound Imaging of Upper Extremity Arterial Aneurysms  
Upper Extremity Arterial Access  
Radial Artery Harvest  
Summary  
References

### 14 Ultrasound Evaluation Before and After Hemodialysis Access

Abstract  
Introduction  
Basic Concepts of Hemodialysis Access  
Descriptive Terminology  
Vascular Mapping Before Hemodialysis Access Placement  
Arteriovenous Fistula Maturity Assessment  
Graft Evaluation  
Acknowledgments  
References

### 15 Ultrasound Assessment of Lower Extremity Arteries

Abstract  
Introduction  
Instrumentation  
Duplex Ultrasound Technique  
Classification of Disease  
Clinical Applications  
Summary  
References

### 16 Ultrasound Assessment During and After Carotid Peripheral Interventions

Introduction  
Intraprocedural Duplex Ultrasound Assessment  
Surveillance After Carotid Intervention  
Surveillance After Peripheral Intervention  
Infrainguinal Bypass Assessment  
Synthetic Grafts  
Duplex-Monitored Peripheral Angioplasty/Stent Placement  
Post-Procedural Follow-Up  
Summary  
References

### 17 Ultrasound in the Assessment and Management of Arterial Emergencies

Abstract  
Introduction  
Ruptured Abdominal Aortic Aneurysm  
Carotid Artery Stenosis/Thrombosis  
Carotid Artery Dissection  
Acute Lower Extremity Ischemia  
Femoral Pseudoaneurysm  
Traumatic Arteriovenous Fistula  
Penetrating Arterial Fistula  
Summary  
References

## Section 4—Extremity Veins

### 18 Extremity Venous Anatomy and Technique for Ultrasound Examination

Introduction  
Anatomy of the Lower Extremity  
Venous Duplex Imaging Examination Technique and Protocol  
The Venous Anatomy and Examination Protocol  
Anatomy of the Upper Extremity

Upper Extremity Examination Protocol  
Characterization of Thrombus  
Summary  
References

## 19 Ultrasound Diagnosis of Lower Extremity Venous Thrombosis

Abstract  
Introduction  
Prevalence, Etiology, and Risk Factors  
Venous Anatomy of the Lower Extremities  
Technique  
Ultrasound Findings  
Pitfalls  
May-Thurner Syndrome  
Alternative Diagnoses/Incidental Findings  
Summary  
References

## 20 Risk Factors and the Role of Ultrasound in the Management of Extremity Venous Disease

Abstract  
Introduction  
Acute Deep Vein Thrombosis Etiology and Risk Factors  
Anticoagulation and Thrombolysis in the Management of Venous Thrombosis  
Acute Deep Vein Thrombosis of Specific Extremity Veins  
Lower Extremity Duplex Scanning in Suspected Pulmonary Embolism  
Lower Extremity Superficial Thrombophlebitis  
Upper Extremity Venous Thrombosis  
Sequelae of Deep Vein Thrombosis  
Primary Varicose Veins  
Preprocedure Venous Mapping  
Summary  
References

## 21 Ultrasound Diagnosis of Venous Insufficiency

Abstract  
Introduction  
Lower Limb Venous Anatomy  
Pathophysiology of Venous Insufficiency  
Duplex Diagnosis of Venous Insufficiency  
Quantitative Measurement of Venous Incompetence  
The Role of Sonography in the Treatment of Chronic Venous Insufficiency  
Summary  
Acknowledgments  
References

## 22 Nonvascular Findings Encountered During Venous Sonography

Introduction  
Soft Tissue Edema  
Lymphedema  
Hematoma  
Muscle Injury  
Lymph Nodes  
Popliteal (Baker) Cysts  
Joint Effusion

Infection  
Soft Tissue Tumors  
References

## Section 5—Abdomen and Pelvis

### 23 Anatomy and Normal Doppler Signatures of Abdominal Vessels

Abstract  
Introduction  
Abdominal Aorta  
Celiac Artery  
Splenic Artery  
Hepatic Artery  
Superior and Inferior Mesenteric Arteries  
Renal Arteries  
Portal Venous System  
Hepatic Veins  
Inferior Vena Cava  
Renal Veins  
References

### 24 Ultrasound Assessment of the Abdominal Aorta

Abstract  
Introduction  
Anatomy  
Normal Sizes  
Normal Doppler Velocity Profiles  
Pathologic States  
Abdominal Aortic Aneurysms  
Ultrasound Examination  
Definition of an Aneurysm: Size Thresholds  
Growth Rate of Aneurysms  
Common Iliac Artery Aneurysms  
Aneurysm Complications  
Varied Pathologies  
Postsurgical Assessment  
Summary  
References

### 25 Ultrasound Assessment Following Endovascular Aortic Aneurysm Repair

Abstract  
Introduction  
Overview of EVAR  
Ultrasound Examination  
Endoleak Detection  
Aneurysm Size  
Endovascular Graft Deformity and Native Artery Complications  
Contrast-Enhanced Ultrasound  
Aneurysm Sac Pressure Measurement  
Surveillance Intervals  
Summary  
References

### 26 Doppler Ultrasound of the Mesenteric Vasculature

Abstract  
Introduction  
Anatomy, Physiology, and Natural History of Bowel Ischemia  
Technique

Examination Protocol  
Diagnostic Criteria  
Roles of Duplex Ultrasound in Surveillance Following  
Mesenteric Revascularization (Stent Angioplasty and  
Bypass Graft Assessment)  
Keys to a Successful Examination  
Pitfalls  
Other Mesenteric Artery Pathologies  
Summary  
Acknowledgments  
References

## 27 Ultrasound Assessment of the Hepatic Vasculature

Abstract  
Introduction  
Technique and Normal Hemodynamics  
Portal Hypertension  
Portal Vein Thrombosis  
Portal Vein Stenosis  
Arteriportal Fistulas  
Portal Vein Gas  
Hepatic Vein Obstruction  
Transjugular Intrahepatic Portosystemic Shunts  
Congenital Intrahepatic Portosystemic Shunts  
Hepatic Artery Stenosis/Occlusion  
Hepatic Artery Pseudoaneurysm  
References

## 28 Duplex Ultrasound of Native Renal Vasculature

Abstract  
Introduction  
Anatomy  
Principles of Examination  
Technique  
Protocol  
Vascular Disorders  
Summary  
Acknowledgments  
References

## 29 Duplex Ultrasound Evaluation of the Uterus and Ovaries

Abstract  
Introduction  
Technical Issues  
Normal Anatomy and Hemodynamics  
Current Applications  
Ectopic Pregnancy  
Pseudogestation  
Value of Doppler for Evaluation of Ectopic Pregnancy  
Pitfalls  
Summary  
Acknowledgments  
References

## 30 Duplex Ultrasound Evaluation of the Male Genitalia

Abstract  
Introduction  
The Scrotum

The Penis  
Summary  
References

## 31 Evaluation of Organ Transplants

Abstract  
Introduction  
Renal Transplantation  
Hepatic Transplantation  
References

## Section 6—Trends in Ultrasound Vascular Imaging

### 32 Credentialing, Accreditation, and Quality in the Vascular Laboratory

Abstract  
Introduction  
Credentialing  
Accreditation  
Quality Improvement  
Appropriateness of Vascular Testing  
Efforts to Standardize Vascular Testing and Improve Quality at  
the National Level  
Summary  
References

### 33 Ultrasound Screening for Vascular Disease

Abstract  
Introduction  
Definition and Types of Screening  
Screening for Asymptomatic Carotid Stenosis  
Abdominal Aneurysm Screening  
Screening for Cardiovascular Disease: Risk and Subclinical  
Cardiovascular Disease  
Summary  
References

### 34 Correlative Imaging

Abstract  
Overview  
Catheter-Based Arteriography  
Magnetic Resonance Angiography  
Computed Tomographic Angiography  
Overview and Correlative Findings  
Summary  
References

### 35 Ultrasound Contrast Agents in Vascular Disease

Abstract  
Introduction  
Technical Aspects of Contrast-Enhanced Ultrasound  
Safety of Ultrasound Contrast Agents  
Clinical Applications of Ultrasound Contrast Agents  
Trauma  
Future Considerations  
Summary  
References

Index

# DAVIES

LIFE IS A TEST...

# Vascular



An Interactive Q&A Review for the ARDMS Specialty Exam

CD-ROM Version 3.1

Donald P. Ridgway, RVT

Cindy A. Owen, RT, RVT, RDMS

Barton A. Bean, RVT

D. E. Strandness, Jr., MD

**DAVIES**

Registry Reviews & Study Aids



**Welcome!**

Thanks for choosing the Davies CD-ROM Mock Exam for Vascular Technology. This powerful learning resource will make a significant difference in your performance on the registry exam.

With this CD you can:

- Study and learn by selecting **Learn Mode**
- Test yourself by selecting **Test Mode**
- Score yourself by viewing a complete **Results Analysis** at the end of your test
- Focus on image-based questions by selecting **Image Library**
- Earn 15 hours of SVU-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- How to use this CD
- 1-2-3 Step Ultrasound Education & Test Prep
- Earn CME credit
- The ARDMS exam content outline
- System requirements
- Legal notices


[Test myself](#)

[Study and learn](#)

[Browse the Image Library](#)

[Earn continuing medical education \(CME\) credit](#)

[Review test-prep and exam-taking strategies](#)

[Access the Contact Directory](#)

**DAVIES**

Registry Reviews &amp; Study Aids

## Welcome!

Thanks for choosing the Davies CD-ROM Mock Exam for Vascular Technology. This powerful learning resource will make a significant difference in your performance on the registry exam.

With this CD you can:

- Study and learn by selecting **Learn Mode**
- Test yourself by selecting **Test Mode**
- Score yourself by viewing a complete **Results Analysis** at the end of your test
- Focus on image-based questions by selecting **Image Library**
- Earn 15 hours of SVU-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- How to use this CD
- 1-2-3 Step Ultrasound Education & Test Prep
- Earn CME credit
- The ARDMS exam content outline
- System requirements
- Legal notices

## Create a Standard ARDMS-Style Test

This program automatically randomizes questions on all ARDMS exam topics in the same ratio as on the real exam.

**Select the number of questions you want . . .**

60  
 120  
 180

Create My Test

Clear My Selection

Return To Start

## Need to focus on specific exam topics?

You can create your own test by selecting as many as three ARDMS exam topics and specifying the number of questions you want in each topic area. To do so, click the purple button below.

This is an excellent way to focus on topics you need to work on.

[Click here to design your own custom test.](#)

## Welcome!

Thanks for choosing the Davies CD-ROM Mock Exam for Vascular Technology. This powerful learning resource will make a significant difference in your performance on the registry exam.

### With this CD you can:

- Study and learn by selecting Learn Mode
- Test yourself by selecting Test Mode
- Score yourself by viewing a complete Results Analysis at the end of your test
- Focus on image-based questions by selecting Image Library
- Earn 15 hours of SVU-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- How to use this CD
- 1-2-3 Step Ultrasound Education & Test Prep
- Earn CME credit
- The ARDMS exam content outline
- System requirements
- Legal notices

## Design My Own Test

Create your own test by selecting as many as three ARDMS exam topics and specifying the number of questions you want in each topic area.

This is an excellent way to focus on topics you need to work on.

### Click below to select one or more exam topics . . .




Total number of questions on the topic you selected.

0 Enter the number of questions you want.

Total number of questions on the topic you selected.

0 Enter the number of questions you want.

Total number of questions on the topic you selected.

0 Enter the number of questions you want.

[Start Page](#)[Learn Mode](#)[Test Mode](#)[Results Analysis](#)[Help](#)[About This CD](#)[Quit](#)

160 The ultrasound image below shows an internal carotid artery with:


- A A calcified plaque
- B An ulcerated lesion
- C A normal arterial wall
- D An intraplaque hemorrhage
- E A homogeneous plaque

[Next](#)


160 of 665

[Prior](#)[Bookmark/Skip](#)[Pause](#)[First](#)[Last](#)[Finish](#)

Question Image 1

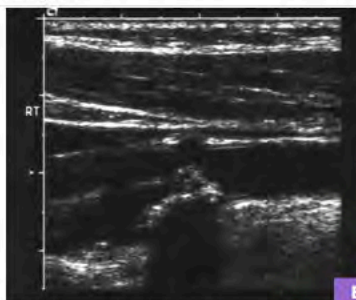[Enlarge](#)

Question Image 2

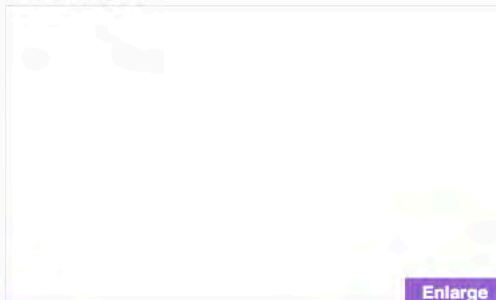[Enlarge](#)

Exam Topic 1: Cerebrovascular (25%–35% of ARDMS exam)

[References](#)[Tutorial](#)[More Images](#)[More Q & A This Topic](#)[Other Learning Tools](#)

Start Page

Learn Mode

Test Mode

Results Analysis

Help

About This CD

Quit

160 The ultrasound image below shows an internal carotid artery with:

- A** A calcified plaque
- B** An ulcerated lesion
- C** A normal arterial wall
- D** An intraplaque hemorrhage
- E** A homogeneous plaque

**RIGHT!**

You chose the correct answer, A: A calcified plaque

*Note the brightly echogenic plaque and acoustic shadowing.*

**References**

Zwiebel WJ: *Introduction to Vascular Ultrasonography*, 5th edition. Philadelphia, Elsevier Saunders, 2005, pp 157–167.

**More Q&A Information**

This question belongs to the following ARDMS exam topic:

I. Cerebrovascular [25–35% of ARDMS exam]

To focus only on this or other specific exam topics:

Next

160 of 665

Prior

Bookmark/Skip


Pause

First

Last


Finish

Explanation Image 1

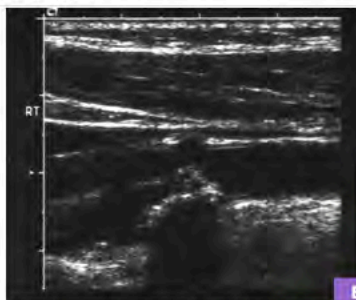

Enlarge

Explanation Image 2

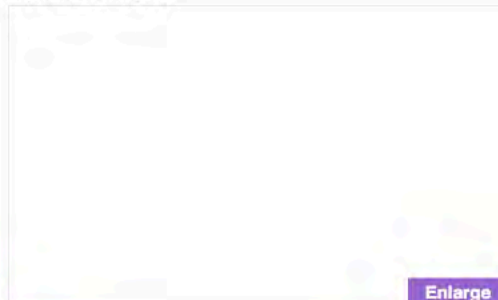

Enlarge

Exam Topic I. Cerebrovascular (25%–35% of ARDMS exam)

References

Tutorial

More Images

More Q &amp; A This Topic

Other Learning Tools

Start Page

Learn Mode

Test Mode

Results Analysis

Help

About This CD

Quit

24 The great saphenous vein:

- A** Originates along the medial dorsum of the foot
- B** Passes superiorly, anterior to the medial malleolus
- C** Is accompanied by the saphenous nerve
- D** Receives tributaries from all surfaces of the lower extremity
- E** All are correct

**Tutorial**

Anatomy of the Lower Extremity Veins

by  
David S. Sumner, MD

Unlike their arterial counterparts, the veins of the leg are divided into two systems: the superficial and the deep (Figures 1 and 2). [Click "More Images" on the purple navigation bar below to view figures.] The two principal veins in the superficial system are the greater (or long) saphenous vein and the lesser (or short) saphenous vein. Beginning just anteriorly and laterally to the medial malleolus, the great saphenous vein courses up the medial aspect of the calf and thigh to enter the common femoral vein at the groin (Figure 3—click "More Images"). The small saphenous vein, which is found posterior to the lateral malleolus, runs up the posterior aspect of the calf and terminates

Next

24 of 665

Prior

Bookmark/Skip

Pause

First

Last

Finish

Explanation Image 1

Explanation Image 2

**Doppler Evaluation of the Lower Extremity Veins**

Click this screen to display QuickTime's video controls, including the play button. The QuickTime controls allow you to start, pause, stop, fast forward, and rewind the tutorial. Click the Play button (▶) to start your video tutorial.

Copyright © 2007 by Davies Publishing. All rights reserved.

Enlarge

Enlarge

Exam Topic I-IV. Anatomy, Physiology &amp; Hemodynamics (4%–18% of ARDMS exam)

References

Tutorial

More Images

More Q &amp; A This Topic

Other Learning Tools

## RESULTS FOR CURRENT TEST

Date: 5/12/2014

Time: 2:45:57 PM

Your results appear below in three categories:

- (1) Overall results.
- (2) Results by exam topic.
- (3) Results by individual questions.

This analysis also tells you how many questions (if any) you failed to answer and for which questions you changed your first answer. From this page you can also print your results, review wrong answers in Learn Mode, repeat the same test again, return to the Start Page for another activity, or Quit.

## OVERALL RESULTS

Total Questions = 60

Total percentage correct = 18.33% / 11 of 60

Target time / actual time used = 60 minutes / 1.02 minutes

Number of unanswered questions = 0

## SUBJECT AREA RESULTS

I-IV. Anatomy, Physiology & Hemodynamics = 33.33% / 1 of 3 (Unanswered 0)

I. Cerebrovascular = 15.38% / 2 of 13 (Unanswered 0)

II. Venous = 15.38% / 2 of 13 (Unanswered 0)

III. Peripheral Arterial = 18.18% / 2 of 11 (Unanswered 0)

[Start Page](#)[Learn Mode](#)[Test Mode](#)[Results Analysis](#)[Help](#)[About This CD](#)[Quit](#)

## SUBJECT AREA RESULTS

I-IV. Anatomy, Physiology & Hemodynamics = 33.33% / 1 of 3 (Unanswered 0)

I. Cerebrovascular = 15.38% / 2 of 13 (Unanswered 0)

II. Venous = 15.38% / 2 of 13 (Unanswered 0)

III. Peripheral Arterial = 18.18% / 2 of 11 (Unanswered 0)

IV. Abdomen and Visceral = 33.33% / 1 of 3 (Unanswered 0)

V. Miscellaneous Conditions & Tests = 33.33% / 1 of 3 (Unanswered 0)

VI. Quality Assurance = 0% / 0 of 2 (Unanswered 0)

VII. Image-Based Questions = 10% / 1 of 10 (Unanswered 0)

VIII. Physiology & Fluid Dynamics = 50% / 1 of 2 (Unanswered 0)

## INDIVIDUAL QUESTION RESULTS

For your convenience, the following question numbers correspond to those in the book version of this mock exam. Use the CD and the book together for best results.

Question ID VTR583 = WRONG VII. Image-Based Questions

Question ID VTR426 = WRONG III. Peripheral Arterial (20%–30% of ARDMS exam)

Question ID VTR570 = WRONG VII. Image-Based Questions

Question ID VTR11 = WRONG I-IV. Anatomy, Physiology & Hemodynamics (4%–18% of ARDMS exam)

Question ID VTR446 = WRONG IV. Abdomen and Visceral (5%–15% of ARDMS exam)

Question ID VTR400 = WRONG III. Peripheral Arterial (20%–30% of ARDMS exam)

Question ID VTR403 = WRONG III. Peripheral Arterial (20%–30% of ARDMS exam)

Question ID VTR502 = WRONG VI. Quality Assurance (3%–5% of ARDMS exam)

Question ID VTR615 = WRONG VII. Image-Based Questions

Question ID VTR531 = WRONG VII. Image-Based Questions

Question ID VTR537 = WRONG VII. Image-Based Questions

Question ID VTR137 = WRONG I. Cerebrovascular (25%–35% of ARDMS exam)

[Review Missed Questions](#)[Review Bookmarked Questions](#)[Review All Questions](#)[Copy Results To Clipboard](#)[Print These Results](#)



**Welcome!**

Thanks for choosing the Davies CD-ROM Mock Exam for Vascular Technology. This powerful learning resource will make a significant difference in your performance on the registry exam.

**With this CD you can:**

- Study and learn by selecting **Learn Mode**
- Test yourself by selecting **Test Mode**
- Score yourself by viewing a complete **Results Analysis** at the end of your test
- Focus on image-based questions by selecting **Image Library**
- Earn 15 hours of SVU-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- How to use this CD
- 1-2-3 Step Ultrasound Education & Test Prep
- Earn CME credit
- The ARDMS exam content outline
- System requirements
- Legal notices

**CME Applicant Information**

Name

Current Credentials

Date of Birth

mm/dd/yyyy

Home Address

City/State/ZIP (include postcode and country, if outside the USA)

Phone

Fax

eMail

ARDMS #

ARRT #

SDMS #

CCI #

**IMPORTANT NOTE:**

Please double check all of your information to ensure smooth, error free processing of your CME application and quiz.

Continue

Clear All Fields

Return To Start

**Copyright © 2019 Davies Publishing Inc.  
All rights reserved.**

**Licensed for individual use only,  
subject to the terms of the License Agreement.  
It is illegal to copy any part of this CD.  
Site licenses are available for ultrasound education.**

**We appreciate your kind cooperation.**

**Questions?  
Call toll-free 1-877-792-0005**

**Visit us online @ [DaviesPublishing.com](http://DaviesPublishing.com)**

# DAVIES

LIFE IS A TEST...

# Physics



An Interactive Q&A Review for the ARDMS SPI Exam

CD-ROM Version 3.1

Cindy A. Owen, RT, RVT, RDMS

James A. Zagzebski, PhD



## Welcome!

Thanks for choosing the Davies CD-ROM Mock Exam for Ultrasound Physics. This powerful learning resource will make a significant difference in your performance on the registry exam.

### With this CD you can:

- Study and learn by selecting Learn Mode
- Test yourself by selecting Test Mode
- Score yourself by viewing a complete Results Analysis at the end of your test
- Focus on image-based questions by selecting Image Library
- Earn 15 hours of SDMS-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- How to use this CD
- 1-2-3 Step Ultrasound Education & Test Prep
- Earn CME credit
- The ARDMS exam content outline
- System requirements
- Legal notices

**Test myself**

**Study and learn**

**Browse the Image Library**

**Earn continuing medical education (CME) credit**

**Review test-prep and exam-taking strategies**

**Access the Contact Directory**

**DAVIES**

Registry Reviews &amp; Study Aids

## Welcome!

Thanks for choosing the Davies CD-ROM Mock Exam for Ultrasound Physics. This powerful learning resource will make a significant difference in your performance on the registry exam.

With this CD you can:

- Study and learn by selecting Learn Mode
- Test yourself by selecting Test Mode
- Score yourself by viewing a complete Results Analysis at the end of your test
- Focus on image-based questions by selecting Image Library
- Earn 15 hours of SDMS-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- How to use this CD
- 1-2-3 Step Ultrasound Education & Test Prep
- Earn CME credit
- The ARDMS exam content outline
- System requirements
- Legal notices

## Create a Standard ARDMS-Style Test

This program automatically randomizes questions on all ARDMS exam topics in the same ratio as on the real exam.

**Select the number of questions you want . . .**

60  
 120  
 180

Create My Test

Clear My Selection

Return To Start

### Need to focus on specific exam topics?

You can create your own test by selecting as many as three ARDMS exam topics and specifying the number of questions you want in each topic area. To do so, click the purple button below.

This is an excellent way to focus on topics you need to work on.

[Click here to design your own custom test.](#)

## Welcome!

Thanks for choosing the Davies CD-ROM Mock Exam for Ultrasound Physics. This powerful learning resource will make a significant difference in your performance on the registry exam.

### With this CD you can:

- Study and learn by selecting Learn Mode
- Test yourself by selecting Test Mode
- Score yourself by viewing a complete Results Analysis at the end of your test
- Focus on image-based questions by selecting Image Library
- Earn 15 hours of SDMS-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- How to use this CD
- 1-2-3 Step Ultrasound Education & Test Prep
- Earn CME credit
- The ARDMS exam content outline
- System requirements
- Legal notices

## Design My Own Test

Create your own test by selecting as many as three ARDMS exam topics and specifying the number of questions you want in each topic area.

This is an excellent way to focus on topics you need to work on.

### Click below to select one or more exam topics . . .




Total number of questions on the topic you selected.

0 Enter the number of questions you want.

Total number of questions on the topic you selected.

0 Enter the number of questions you want.

Total number of questions on the topic you selected.

0 Enter the number of questions you want.

[Create My Test](#)
[Clear My Selections](#)
[Return To Start](#)

Start Page

Learn Mode

Test Mode

Results Analysis

Help

About This CD

Quit

48 Which of the following sound interactions produced the shadows indicated by the arrows on this image?

- A Reflection
- B Diffraction
- C Diffusion
- D Rayleigh scattering
- E Refraction

Next

48 of 600

Prior

Bookmark/Skip

Pause

First

Last


Finish

Question Image 1

Question Image 2



Enlarge

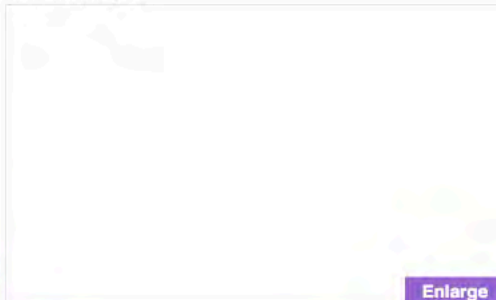

Enlarge

Exam Topic II. Physical Principles (20% of ARDMS exam)

References

Tutorial

More Images

More Q &amp; A This Topic

Other Learning Tools

Start Page

Learn Mode

Test Mode

Results Analysis

Help

About This CD

Quit

48 Which of the following sound interactions produced the shadows indicated by the arrows on this image?

A Reflection

B Diffraction

C Diffusion

D Rayleigh scattering

E Refraction

**RIGHT!**

You chose the correct answer, E:  
Refraction

*The **edge shadowing** seen with curved interfaces is caused by refraction of the sound beam. This is the same physical principle that can be observed by placing a pencil in a partially filled glass of water. When viewed from the side, the pencil appears bent.*

Next

48 of 600

Prior

Bookmark/Skip

Pause

First

Last


Finish

Explanation Image 1

Explanation Image 2



Enlarge

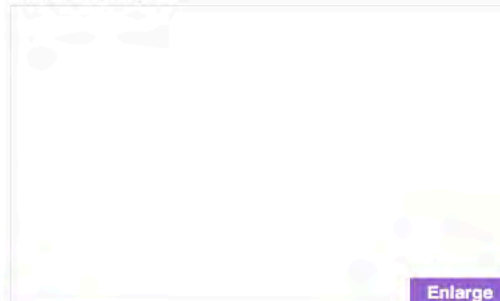

Enlarge

### More Q&A Information

This question belongs to the following ARDMS exam topic:

II. Physical Principles (20% of ARDMS exam)

To focus only on this or other specific exam topics:

Exam Topic II. Physical Principles (20% of ARDMS exam)

References

Tutorial

More Images

More Q &amp; A This Topic

Other Learning Tools



## RESULTS FOR CURRENT TEST

Date: 5/12/2014

Time: 3:40:07 PM

Your results appear below in three categories:

- (1) Overall results.
- (2) Results by exam topic.
- (3) Results by individual questions.

This analysis also tells you how many questions (if any) you failed to answer and for which questions you changed your first answer. From this page you can also print your results, review wrong answers in Learn Mode, repeat the same test again, return to the Start Page for another activity, or Quit.

## OVERALL RESULTS

Total Questions = 60

Total percentage correct = 13.33% / 8 of 60

Target time / actual time used = 60 minutes / 1.41 minutes

Number of unanswered questions = 0

## SUBJECT AREA RESULTS

- I. Patient Care, Safety, and Communication = 0% / 0 of 3 (Unanswered 0)
- II. Physical Principles = 16.67% / 2 of 12 (Unanswered 0)
- III. Ultrasound Transducers = 9.09% / 1 of 11 (Unanswered 0)
- IV. Pulse-Echo Instrumentation = 16.67% / 3 of 18 (Unanswered 0)

**SUBJECT AREA RESULTS**

- I. Patient Care, Safety, and Communication = 0% / 0 of 3 (Unanswered 0)
- II. Physical Principles = 16.67% / 2 of 12 (Unanswered 0)
- III. Ultrasound Transducers = 9.09% / 1 of 11 (Unanswered 0)
- IV. Pulse-Echo Instrumentation = 16.67% / 3 of 18 (Unanswered 0)
- V. Doppler Instrumentation and Hemodynamics = 8.33% / 1 of 12 (Unanswered 0)
- VI. Quality Assurance/Quality Control of Equipment = 25% / 1 of 4 (Unanswered 0)

**INDIVIDUAL QUESTION RESULTS**

For your convenience, the following question numbers correspond to those in the book version of this mock exam. Use the CD and the book together for best results.

- Question ID SPI211 = WRONG III. Ultrasound Transducers (20% of ARDMS exam)
- Question ID SPI87 = WRONG II. Physical Principles (20% of ARDMS exam)
- Question ID SPI113 = WRONG II. Physical Principles (20% of ARDMS exam)
- Question ID SPI110 = WRONG II. Physical Principles (20% of ARDMS exam)
- Question ID SPI588 = RIGHT VI. Quality Assurance/Quality Control of Equipment (5% of ARDMS exam)
- Question ID SPI250 = WRONG III. Ultrasound Transducers (20% of ARDMS exam)
- Question ID SPI121 = WRONG II. Physical Principles (20% of ARDMS exam)
- Question ID SPI242 = RIGHT III. Ultrasound Transducers (20% of ARDMS exam)
- Question ID SPI590 = WRONG VI. Quality Assurance/Quality Control of Equipment (5% of ARDMS exam)
- Question ID SPI463 = WRONG V. Doppler Instrumentation and Hemodynamics (20% of ARDMS exam)
- Question ID SPI507 = WRONG V. Doppler Instrumentation and Hemodynamics (20% of ARDMS exam)
- Question ID SPI407 = WRONG IV. Pulse-Echo Instrumentation (30% of ARDMS exam)
- Question ID SPI553 = RIGHT V. Doppler Instrumentation and Hemodynamics (20% of ARDMS exam)
- Question ID SPI539 = WRONG V. Doppler Instrumentation and Hemodynamics (20% of ARDMS exam)
- Question ID SPI348 = WRONG IV. Pulse-Echo Instrumentation (30% of ARDMS exam)

**Welcome!**

Thanks for choosing the Davies CD-ROM Mock Exam for Ultrasound Physics. This powerful learning resource will make a significant difference in your performance on the registry exam.

**With this CD you can:**

- Study and learn by selecting **Learn Mode**
- Test yourself by selecting **Test Mode**
- Score yourself by viewing a complete **Results Analysis** at the end of your test
- Focus on image-based questions by selecting **Image Library**
- Earn 15 hours of SDMS-approved CME credit

To start immediately, choose any of the activities on the right. To learn more about this powerful learning tool from Davies, click any of the following topics:

- **How to use this CD**
- **1-2-3 Step Ultrasound Education & Test Prep**
- **Earn CME credit**
- **The ARDMS exam content outline**
- **System requirements**
- **Legal notices**

**CME Applicant Information**

Name

Current Credentials

Date of Birth

mm/dd/yyyy

Home Address

City/State/ZIP (Include postcode and country if outside the USA)

Phone

Fax

eMail

ARDMS #

ARRT #

SDMS #

CCI #

**IMPORTANT NOTE:**

Please double check all of your information to ensure smooth, error free processing of your CME application and quiz.

Continue

Clear All Fields

Return To Start

**Copyright © 2019 Davies Publishing Inc.  
All rights reserved.**

**Licensed for individual use only,  
subject to the terms of the License Agreement.  
It is illegal to copy any part of this CD.  
Site licenses are available for ultrasound education.**

**We appreciate your kind cooperation.**

**Questions?  
Call toll-free 1-877-792-0005**

**Visit us online @ [DaviesPublishing.com](http://DaviesPublishing.com)**